

# ULT Freezer Solutions



● Headquarters ● Global warehouse ● Training center ● Global subsidiary ● Manufacturing base



# CONTENT

**-150°C Cryo Freezer** ..... 01

**-86°C ULT Freezer** ..... 05

- TwinCool Frequency Conversion ULT Freezer ..... 05
- TwinCool ULT Freezer ..... 11
- Frequency Conversion ULT Freezer ..... 19
- Water-Cooled ULT Freezer ..... 27
- Standard Low Energy ULT Freezer with LCD Touchscreen --- 29
- Standard Low Energy ULT Freezer with LED Display ..... 33
- Under Counter ULT Freezer ..... 40
- Controlled Rate Freezer ..... 49
- Backup System ..... 51

**Consumables Solutions for ULT Freezer** ..... 53

- Freezer Storage Racks ..... 55
- Barcode Cryogenic Tubes and Storage Boxes ..... 57
- Universal Cryogenic Tubes and Storage Boxes ..... 58
- High Speed Batch Scanner ..... 59
- Single Tube Scanner ..... 61
- Batch Cover Opener ..... 62
- Low Temperature Label Printer ..... 63

**Product Portfolio** ..... 65





**-150°C  
Cryo Freezer**



A mechanical cryogenic freezer is used for the long-term cryopreservation of samples such as cell cultures, tissue samples, bacteria, viruses and non-biological products such as electronics. Applications include cord blood preservation, biobanks, reproductive procedures, blood banks, life science research, biomedical engineering and fisheries. Mechanical cryopreservation is a safe alternative to LN<sub>2</sub> and can reduce or eliminate the reliance on LN<sub>2</sub>, reducing safety concerns.



#### Energy Efficient Hydrocarbon Refrigeration

- Includes specialist industrial compressor, hydrocarbon refrigerant, dual-compressors 4-stage cascade refrigeration system, which delivers excellent energy-efficiencies
- Energy-efficient fan design with intelligent start/stop



#### Environmentally Friendly

- **100%** HC refrigerant truly environmental-friendly refrigeration system with zero global warming potential (GWP)
- LBA polyurethane foam has low GWP and is significantly lower than most HFCs commonly used, offering better environmental protection



#### Superior Uniformity

The temperature **uniformity is ±5°C** with the setting temperature of -150°C



#### Low Noise

Noise-reduction features and noise reduction design lowers to noise to 55 dB (A)



#### Safe, Secure and Reliable

Multiple alarms, including high and low temperature, power failure, door ajar, hot condenser, low battery, sensor error, and high ambient temperature. The standard ambient temperature ranges from 10°C to 32°C. The pressure state of the system can be automatically adjusted by the pressure controller to prevent the pressure from being too high to effect the freezer operation.



## Specifications


Model		DW-150W209	DW-150W209A	
Technical Data	Cabinet Type	Chest	Chest	
	Climate Class	N	N	
	Cooling Type	Direct cooling	Direct cooling	
	Defrost Mode	Manual	Manual	
	Refrigerant	HC	HC	
	Sound Level (dB(A))	55	55	
Performance	Cooling Performance (°C)	-150	-150	
	Temperature Range (°C)	-140--150	-140--150	
Control	Controller	Microprocessor	Microprocessor	
	Display	LCD	LCD	
Electrical Data	Power Supply (V/Hz)	220-240/50	208-230/60	
	Electrical Current (A)	10	12	
Dimensions	Capacity (L/Cu.Ft)	209/7.4	209/7.4	
	Net/Gross Weight	kg	260/295	260/295
		lbs	573.2/650.4	573.2/650.4
	Interior Dimension (W*D*H)	mm	667*462*677	667*462*677
		in	26.3*18.2*26.7	26.3*18.2*26.7
	Exterior Dimension (W*D*H)	mm	1650*860*1010	1650*860*1010
		in	65.0*33.9*39.8	65.0*33.9*39.8
	Packing Dimension (W*D*H)	mm	1810*960*1170	1810*960*1170
		in	71.3*37.8*46.1	71.3*37.8*46.1
	Container Load (20'/40'/40'H)		12/24/24	12/24/24
Functions	High/Low Temperature	Y	Y	
	Hot Condenser	Y	Y	
	Power Failure	Y	Y	
	High/Low Voltage	Y	Y	
	Sensor Error	Y	Y	
	Low Battery	Y	Y	
	High Ambient Temperature	Y	Y	
	Door Ajar	Y	Y	
	Start-up Delay	Y	Y	
	Caster	Y	Y	
	Foot	Y	Y	
	Porthole	Y/1	Y/1	
	Shelves/Inner Doors	-/2	-/2	
Accessories	USB Interface	Y	Y	
	Remote Alarm	Y	Y	
	5V Power Supply Port	Y	Y	
	Temperature Chart Recorder	Optional	Optional	
	RS232/485 Port	NA/Standard	NA/Standard	
	CO <sub>2</sub> Backup System	N/A	N/A	
	LN <sub>2</sub> Backup System (Built-in/External)	Optional/Optional	Optional/Optional	
	Freezing Rack	433026, Quantity: 12 433028, Quantity: 12	433026, Quantity: 12 433028, Quantity: 12	
	Other	Certification	CE	UL

Product appearance and specifications are subject to change without notice



**TwinCool  
Frequency Conversion  
ULT Freezer**

### Product Advantages

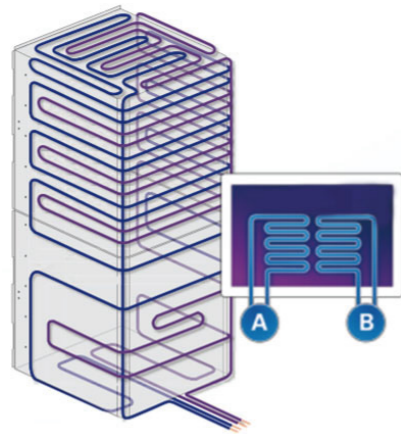

**Dual refrigeration system for optimal reliability and security**  
Two independent hydrocarbon (HC) refrigeration systems, each can maintain -80°C ensuring the safety of stored samples

**World-leading Energy Saving Refrigeration Technology**  
Intelligent frequency conversion technology and HC refrigeration system delivers a 50% reduction in energy usage when compared to traditional HFC refrigeration systems

**Environmental Protection**  
Environmentally friendly hydrocarbon refrigerant and LBA foam material

**Optional IoT System**  
The optional IoT monitors the operating status of the equipment in real time. Equipped with multiple alarm functions and a self-diagnostics system to identify and warn users to ensure sample safety

### Ergonomic Design

**Low Noise Design**  
Optimised systems and noise-reducing cabinet design lower noise out to 43dB for a quieter working experience

**Adjustable Sample Loading Tray (Optional)**  
The adjustable height sample load tray provides users with a convenient place to hold samples while opening and accessing storage compartments


**Large loading capacity**  
Maximum Capacity: 600 boxes (2 inch boxes with 10\*10 configuration) 60,000 samples can be stored

**Intelligent Control**  
10.1-inch high-performance touch screen, sensitive touch operation. Using the optional IoT module, users can check the real-time operating status from an App

**Ergonomic design**  
Document holder for easy storage of documents



**High quality door opening experience**  
The ergonomic chamfer design of the handle ensures smooth operation and maximum durability

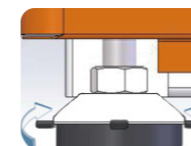

**Adjustable feet**  
For easy positioning and levelling



**Magnetic filter screen**  
Filter screen is easy to clean



DW-86L828BPST

**Microcomputer control**  
Microcontroller with 10.1 inch LCD touch screen displays temperature and other important data. Temperature display precision is 0.1°C

**Multiple alarms**  
Including high/ low temperature alarm, sensor error alarm, hot condensor alarm, power failure alarm, low battery alarm, door ajar alarm, high ambient temperature alarm, etc




**Multiple safety lock**  
Standard mechanical lock, optional electromagnetic lock. Optional fingerprint lock and NFC for optimal security and access




**Networking function**  
Equipped with RS485 interface to connect multiple freezers

### Temperature Curve


Temperature Cool Down Curve

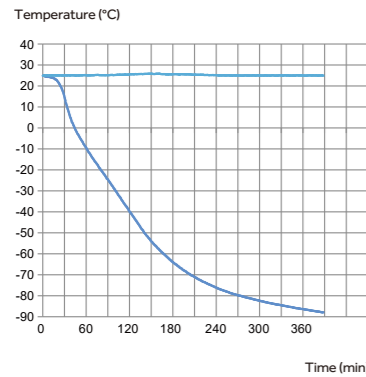

Stability Curve

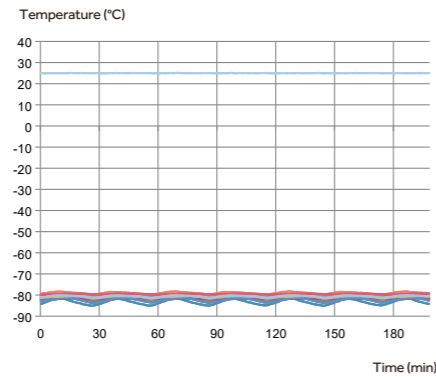

Temperature Recovery Curve

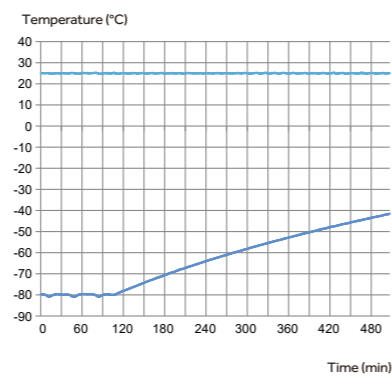

— 3<sup>rd</sup> centre  
— amb temp

— ave temp  
— 1<sup>st</sup> center temp  
— 2<sup>nd</sup> left outer temp  
— 2<sup>nd</sup> right inner temp  
— 3<sup>rd</sup> center temp

— 1<sup>st</sup> left inner temp  
— 1<sup>st</sup> right outer temp  
— 2<sup>nd</sup> center temp  
— 3<sup>rd</sup> left inner temp  
— 3<sup>rd</sup> right outer temp

— 3<sup>rd</sup> centre  
— amb temp

### Product Dimensions



W0	W1	W2	W3	H0	H1	H2	H3	H4	H5	D0	D1	D2
1145mm	1050mm	870mm	854mm	1980mm	1310mm	298mm	303mm	303mm	301mm	998mm	695mm	676mm

### Specifications

Model	DW-86L578BPST	DW-86L728BPST	DW-86L828BPST	
<b>Technical Data</b>				
Cabinet Type	Upright	Upright	Upright	
Climate Class	N	N	N	
Cooling Type	Direct cooling	Direct cooling	Direct cooling	
Defrost Mode	Manual	Manual	Manual	
Refrigerant	HC	HC	HC	
Sound Level (dB(A))	42	42	42	
<b>Performance</b>				
Cooling Performance (°C)	-86	-86	-86	
Temperature Range (°C)	-40--86	-40--86	-40--86	
<b>Control</b>				
Controller	Microprocessor	Microprocessor	Microprocessor	
Display	LCD	LCD	LCD	
<b>Electrical Data</b>				
Power Supply (V/Hz)	100-230/50/60Hz	100-230/50/60Hz	100-230/50/60Hz	
Electrical Current (A)	10.4	10.4	13	
Capacity (L/Cu.Ft)	578/20.4	728/25.7	828/29.2	
Net/Gross Weight (approx)	kg	325/355	350/385	380/410
	lbs	716.5/782.6	771.6/848.8	837.7/914.9
Interior Dimension (W*D*H)	mm	620*716*1310	766*716*1310	870*716*1310
	in	24.4*28.2*51.6	30.2*28.2*51.6	34.3*28.2*51.6
Exterior Dimension (W*D*H)	mm	895*998*1980	1046*998*1980	1145*998*1980
	in	35.2*39.3*78.0	41.2*39.3*78.0	45.1*39.3*78.0
Packing Dimension (W*D*H)	mm	975*1095*2130	1120*1090*2130	1240*1080*2130
	in	38.4*43.1*83.9	44.1*42.9*83.9	48.8*42.5*83.9
Cabinet Width (without handle and hinge)	mm	800	946	1050
Container Load (20'/40'/40'H)		12/24/24	10/22/22	5/11/11
<b>Functions</b>				
High/Low Temperature	Y	Y	Y	
Hot Condenser	Y	Y	Y	
Power Failure	Y	Y	Y	
High/Low Voltage	Y	Y	Y	
Sensor Error	Y	Y	Y	
Low Battery	Y	Y	Y	
High Ambient Temperature	Y	Y	Y	
Door Ajar	Y	Y	Y	
Start-up Delay	Y	Y	Y	
<b>Accessories</b>				
Caster	Y	Y	Y	
Foot	Y	Y	Y	
Porthole	Y/2	Y/2	Y/2	
Shelves/ Inner Doors	3/4	3/4	3/4	
USB Interface	Y	Y	Y	
Remote Alarm	Y	Y	Y	
5V Power Supply Port	Y	Y	Y	
Temperature Chart Recorder	Optional	Optional	Optional	
RS232/485 Port	NA/Standard	NA/Standard	NA/Standard	
CO <sub>2</sub> Backup System	Optional	Optional	Optional	
LN <sub>2</sub> Backup System	Optional	Optional	Optional	
<b>Certifications</b>				
CE	Y	Y	Y	
UL	Y	Y	Y	

# TwinCool ULT Freezer



## Scope of Application

The TwinCool ULT freezer can be used for the storage and protection of valuable samples which require strict and continuous storage conditions, designed to operate even in the event of a compressor failure. Suitable for viruses, pathogens, blood cells and other biological sample cold storage applications found within hospitals, disease control, research institutions and biomedical engineering. Also used to store special materials and other products within electronics and chemicals industries.

## Advanced Hardware System



### Smart Full-size Touch Screen

10-inch touch screen with state-of-art user interface design, coupled with sample management system provides optimal user experience.

## Optional IOT & Sample Management System (Professional Edition)



### Simplified Sample Management Experience

Barcode scanner for simple, effortless and precise identification. Input and retrieve your samples with higher precision and efficiency



### Wireless Monitoring Connectivity

Check the real-time operating status via mobile phones or palmtop, simple and reliable

## Friendly Design



### Safe and secure

Equipped with key lock, padlock and optional electromagnetic lock, with optional fingerprint lock, providing multiple safeguards for sample safety



### Cloud data storage available

Store hundreds of millions of scientific research and sample information in the cloud server

### A Dual independent refrigeration systems for maximum sample safety

The dual refrigeration systems run independently and alternately, both reaching -80°C, such that if one system fails the other will maintain temperature to ensure sample storage safety.

### B High speed refrigeration system for faster pull down and temperature recovery after door opening

Utilizing auto-cascade hydrocarbon (HC) refrigeration technology to deliver fast temperature pull down. From an ambient of 25°C it takes just 180 minutes to reach -80°C. It provides a quick temperature recovery after door opening, with the internal temperatures returning to -75°C within 1 minute, guaranteeing safety of samples.

### C World-leading energy saving refrigeration technology

The HC refrigeration technology coupled with superinsulation, which increases the insulation efficiency by 30%, and a cabinet designed to reduce heat loss ensures an energy efficient freezer. The 578L model has a power consumption of 11 kWh/day and is certified by The National Quality Certification Center for Energy Saving and Environmental Protection.



### Low noise design, reducing the noise down to 53dB

Special noise-reduction design plus super silent compressor technology and energy-saving fan, considerably lowers noise level



### Optimized insulation

Double foaming for both inner and outer doors and five-layer sealing design and optimized super-thick VIP thermal insulation technology, extends temperature holdover time during power failure and increases insulation efficiency by 30%

## Specifications

Model	DW-86L578ST		DW-86L578SAT		DW-86L728ST		DW-86L728SAT			
Technical Data	Cabinet Type	Upright		Upright		Upright		Upright		
	Climate Class	N		N		N		N		
	Cooling Type	Direct cooling		Direct cooling		Direct cooling		Direct cooling		
	Defrost Mode	Manual		Manual		Manual		Manual		
	Refrigerant	HC		HC		HC		HC		
	Sound Level (dB(A))	53	52	53	50	53	53			
Performance	Cooling Performance (°C)	-86		-86		-86		-86		
	Temperature Range (°C)	-40~-86		-40~-86		-40~-86		-40~-86		
Control	Controller	Microprocessor		Microprocessor		Microprocessor		Microprocessor		
	Display	LCD Touchscreen		LCD Touchscreen		LCD Touchscreen		LCD Touchscreen		
Electrical Data	Power Supply (V/Hz)	220-240/50		120/60		208~230/60		208-230/60		
	Electrical Current (A)	10	18	10	10	18	10			
Construction	Capacity (L/Cu.Ft)	578/20.4		578/20.4		728/25.7		728/25.7		
	Net/Gross Weight (approx)	kg	325/355		325/355		350/385		350/385	
		lbs	716.5/782.6		716.5/782.6		771.6/848.8		771.6/848.8	
	Interior Dimension (W*D*H)	mm	620*716*1310		620*716*1310		766*716*1310		766*716*1310	
		in	24.4*28.2*51.6		24.4*28.2*51.6		30.2*28.2*51.6		30.2*28.2*51.6	
	Exterior Dimension (W*D*H)	mm	895*998*1980		895*998*1980		1046*998*1980		1046*998*1980	
		in	35.4*39.3*78.0		35.4*39.3*78.0		41.2*39.3*78.0		41.2*39.3*78.0	
Packing Dimension (W*D*H)	mm	975*1095*2130		975*1095*2130		1120*1090*2130		1120*1090*2130		
	in	38.4*43.1*83.9		38.4*43.1*83.9		44.1*42.9*83.9		44.1*42.9*83.9		
Loading Quantities	Container Load (20'/40'/40'H)	12/24/24		12/24/24		10/22/22		10/22/22		
	High/Low Temperature	Y		Y		Y		Y		
Alarms	Hot Condenser	Y		Y		Y		Y		
	Power Failure	Y		Y		Y		Y		
	High/Low Voltage	Y		Y		Y		Y		
	Sensor Error	Y		Y		Y		Y		
	Low Battery	Y		Y		Y		Y		
	High Ambient Temperature	Y		Y		Y		Y		
	Door Ajar	Y		Y		Y		Y		
	Caster	Y		Y		Y		Y		
	Foot	Y		Y		Y		Y		
	Porthole	Y/2		Y/2		Y/2		Y/2		
Accessories	Shelves/Inner Doors	3/4		3/4		3/4		3/4		
	USB Interface	Y		Y		Y		Y		
	Remote Alarm (Dry contacts)	Y		Y		Y		Y		
	5V Power Supply Port	Y		Y		Y		Y		
	Temperature Recorder	Optional		Optional		Optional		Optional		
	RS485 Port	Y		Y		Y		Y		
	CO <sub>2</sub> Backup System	Optional		Optional		Optional		Optional		
	LN <sub>2</sub> Backup System	Optional		Optional		Optional		Optional		
	Certifications	CE	Y	/	/	Y	/	/		
		UL	/	Y	Y	/	Y	Y		
ENERGY STAR		Y	Y	Y	Y	Y	Y			

Product appearance and specifications are subject to change without notice

The TwinCool ultra-low temperature freezer of Haier Biomedical provides the highest level of protection for your valuable samples.



DW-86L578S

### Safer by Design

Haier Biomedical's ultra low temperature freezers with intelligent TwinCool technology are designed to provide optimal cabinet reliability, longevity, efficiency and sample protection. This super efficient technology also improves the energy efficiency of our third generation ULT freezers and leads the way in terms of product innovation.

### Intelligent TwinCool Refrigeration System

Two independent refrigeration systems are designed to ensure optimal reliability, longevity and efficiency. Depending on the load demands and ambient conditions, one or both refrigeration systems will operate on demand, ensuring the samples are fully protected under the worst possible condition.




#### Maximum Sample Security

TwinCool system means extra insurance for temperature. Each independent refrigeration system can maintain  $-80^{\circ}\text{C}$  separately.



#### Fast Cabinet Pull Down

Fast and efficient cabinet pulldown, it usually takes an average of three hours to reach  $-80^{\circ}\text{C}$  at  $25^{\circ}\text{C}$  ambient. This means the temperature recovery after door opening is excellent ensuring the stored samples are not exposed to undesirable temperatures.

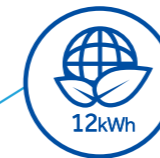

#### Maximum Energy Efficiency

The TwinCool ultra-low temperature system operates with 12 kWh/day.



#### World-leading Energy Saving Refrigeration Technology

The Haier hydrocarbon refrigeration technology uses less than 50% energy compared to traditional CFC refrigerants to reduce the operating cost. The refrigerants do not contain fluorine or chlorine giving it a GWP value of just three, which is better for the environment.



#### Reduced Running Costs

VIP thermal insulation system is designed to significantly reduce heat gain and operating cost.



**Energy Saving Refrigeration**  
High efficiency cooling fans and compressors, combined with hydrocarbon refrigerants, ensure energy savings and long-term sample security



**Multilayered Sealing Structure**  
Multilayered gaskets decrease heat loss and guarantee excellent warm up times in the event of a power failure



**Pressure Equalisation Port**  
Heated Pressure Equalization Port allows users to re-open the main door quickly when entering. Adopts chromium plating, rust-proof



**Improved Handle Design**  
Lockable handle safeguards your precious samples. A padlock can also be added for extra sample safety



**USB Interface**  
Enables users to download historical temperature data for compliance/auditing purposes




**Multi-level Alarms**  
Alarm functions include high, low temperature, sensor error, power failure, high ambient, clean filter and door ajar



DW-86L578S


**DW-86L578S | TYPICAL PERFORMANCE CHARACTERISTICS AT 25°C AMBIENT**

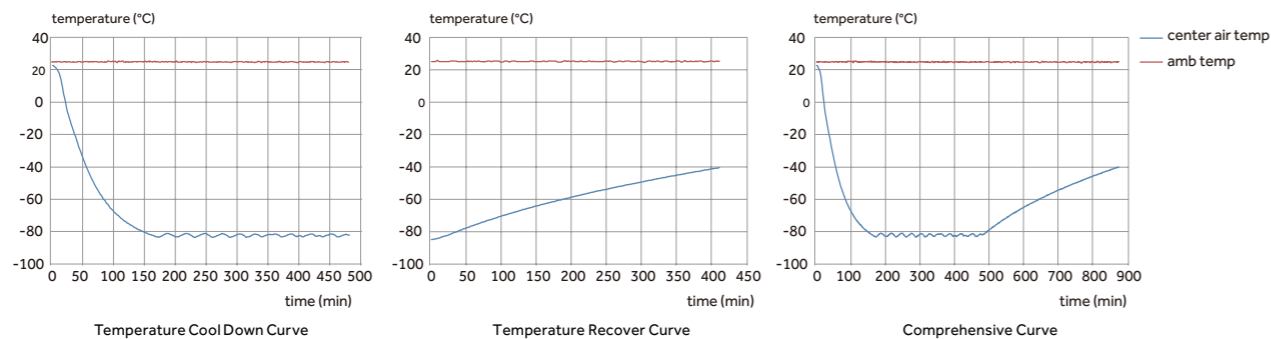

**Specifications**

	Model	DW-86L578S	DW-86L728S	
Technical Data	Cabinet Type	Upright	Upright	
	Climate Class	N	N	
	Cooling Type	Direct Cooling	Direct Cooling	
	Defrost Mode	Manual	Manual	
	Refrigerant	HC	HC	
	Sound Level (dB(A))	53	50	
Performance	Cooling Performance (°C)	-86	-86	
	Temperature Range (°C)	-40~-86	-40~-86	
Control	Controller	Microprocessor	Microprocessor	
	Display	LED	LED	
Electrical Data	Power Supply (V/Hz)	220~240/50	220~240/50	
	Electrical Current (A)	10	10	
Construction	Capacity (L/Cu.Ft)	578/20.4	728/25.7	
	Net/Gross Weight (approx)	kg	325/355	350/385
		lbs	716.5/782.6	771.6/848.8
	Interior Dimension (W*D*H)	mm	620*716*1310	766*716*1310
		in	24.4*28.2*51.6	30.2*28.2*51.6
	Exterior Dimension (W*D*H)	mm	895*998*1980	1046*998*1980
		in	35.4*39.3*78.0	41.2*39.3*78.0
	Packing Dimension (W*D*H)	mm	975*1095*2130	1120*1090*2130
		in	38.4*43.1*83.9	44.1*42.9*83.9
	Loading Quantities	Container load (20'/40'/40'H)	12/24/24	10/22/22
Alarms	High/Low Temperature	Y	Y	
	Hot Condenser	Y	Y	
	Power Failure	Y	Y	
	Sensor Error	Y	Y	
	Low Battery	Y	Y	
	High Ambient Temperature	Y	Y	
	Door Ajar	Y	Y	
Accessories	Caster	Y	Y	
	Foot	Y	Y	
	Porthole	Y/2	Y/2	
	Shelves/Inner doors	3/4	3/4	
	USB Interface	Y	Y	
	Remote Alarm (Dry contacts)	Y	Y	
	5V Power Supply Port	Y	Y	
	Temperature Recorder	Optional	Optional	
	RS485	Y	Y	
	CO <sub>2</sub> Backup System	Optional	Optional	
LN <sub>2</sub> Backup System	Optional	Optional		
Certifications	CE	Y	Y	
	ENERGY STAR	Y	Y	

S suffix - Dual independent refrigeration systems  
Product appearance and specifications are subject to change without notice

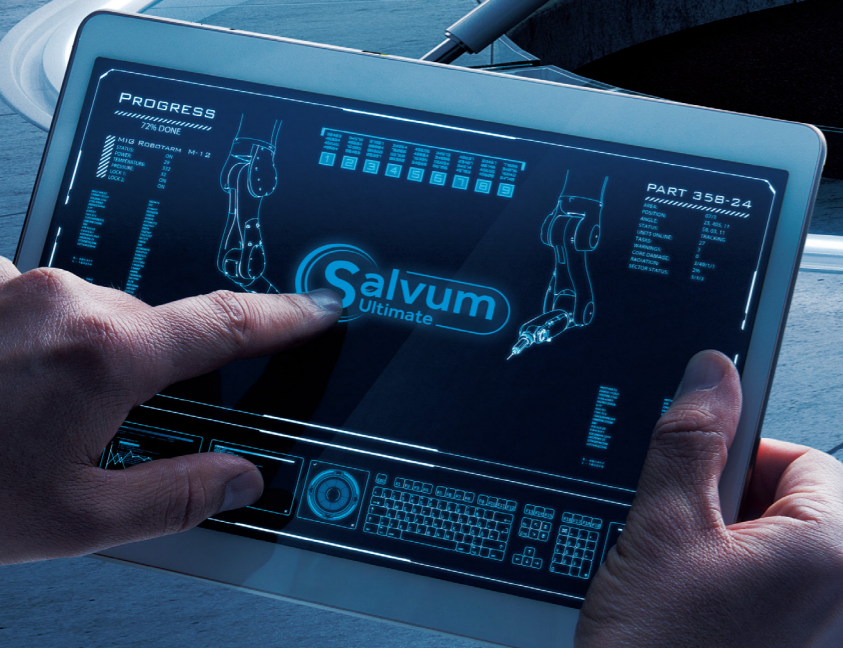

Frequency  
Conversion ULT  
Freezer

## Scope of Application

Applicable for products and samples which require strict storage conditions such as viruses, pathogens, red blood cells, white blood cells, skin, bones, bacteria, semen, biological products, electronics and special materials. Suitable for long term storage applications and compliant with typical storage requirements found in hospitals, disease control and prevention centers, scientific research institutions, biomedical engineering institutes, agriculture/ fishery companies as well as the electronics and chemical industry.

## Advanced Hardware System



### Smart Full-size Touch Screen

10-inch capacitive touch screen, state-of-art UI design coupled with sample management system to provide the optimal user experience.



### HC & Variable Frequency Drive Refrigeration System for Additional Energy Saving

Advanced innovative design delivers excellent energy savings. The energy consumption is down to a single digit.

## Optional IoT Software System



### Simplified Sample Management Experience

Optional barcode scanner for simple, effortless and precise identification. Input and retrieve your samples with higher precision and efficiency



### Wireless Monitoring Connectivity

Check the real-time operating status via mobile phones or palmtop, simple and reliable

### A Quicker product access, identification and retrieval

Instead of a manual system, the one-gun, one-code and one-key operation plus touch screen synchronization means you can access and retrieve your samples within seconds

### B 24-hour sample protection

Using Haier's app and IoT technology the unit can be monitored and can self-diagnose faults, ensuring you are always aware of your unit's status and able to make real-time and informed choices to protect your samples

### C Cloud data storage available

Store hundreds of millions of scientific research and sample information in the cloud server

## Friendly Design



### Safe and secure

Standard equipped with key lock, padlock and electromagnetic lock with optional fingerprint lock, providing multiple safeguards for sample safety



### Low noise design, reducing the noise down to 43.5dB

Special noise-reduction design plus super silent compressor technology and energy-saving fan, considerably lowers noise level



### Super energy efficient with three environmentally friendly innovations

Ultra-low power consumption, down to less than 10kWh/day, ensuring a world leading energy saving performance



### Optimized insulation

Double foaming for both inner and outer doors and five-layer sealing design and optimized super-thick VIP thermal insulation technology, extends temperature holdover time during power failure and increases insulation efficiency by 20%

## Specifications

Model	DW-86L579BPT	DW-86L729BPT	DW-86L829BPT	DW-86L959BPT		
Technical Data	Cabinet Type	Upright	Upright	Upright	Upright	
	Climate Class	N	N	N	N	
	Cooling Type	Direct Cooling	Direct Cooling	Direct Cooling	Direct Cooling	
	Defrost Mode	Manual	Manual	Manual	Manual	
	Refrigerant	HC	HC	HC	HC	
	Sound Level (dB(A))	43.5	43.5	43.5	47	
Performance	Cooling Performance (°C)	-86	-86	-86	-86	
	Temperature Range (°C)	-40~-86	-40~-86	-40~-86	-40~-86	
Control	Controller	Microprocessor	Microprocessor	Microprocessor	Microprocessor	
	Display	Touch Screen LCD	Touch Screen LCD	Touch Screen LCD	Touch Screen LCD	
Electrical Data	Power Supply (V/Hz)	100-230/50/60	100-230/50/60	100-230/50/60 208-230/50/60	100-230/50/60 208-230/50/60	
	Electrical Current (A)	14	14	15 6	20 7	
Dimensions	Capacity (L/Cu.Ft)	579/20.4	729/25.7	829/29.2	959/33.9	
	Net/Gross Weight (approx)	kg	325/355	350/385	380/415	450/485
		lbs	716.5/782.6	771.6/848.8	837.7/914.9	992.1/1069.2
	Interior Dimension (W*D*H)	mm	620*716*1310	766*716*1310	870*716*1310	1016*716*1310
		in	24.4*28.2*51.6	30.2*28.2*51.6	34.3*28.2*51.6	40.0*28.2*51.6
	Exterior Dimension (W*D*H)	mm	895*998*1980	1046*998*1980	1145*998*1980	1296*998*1980
		in	35.2*39.3*78.0	41.2*39.3*78.0	45.1*39.3*78.0	51.0*39.3*78.0
	Packing Dimension (W*D*H)	mm	975*1095*2130	1120*1090*2130	1240*1080*2130	1410*1160*2185
		in	38.4*43.1*83.9	44.1*42.9*83.9	48.8*42.5*83.9	55.5*45.7*86.0
	Container load (20'/40'/40'H)	12/24/24	10/22/22	5/11/11	5/10/10	
Functions	High/Low Temperature	Y	Y	Y	Y	
	Hot Condenser	Y	Y	Y	Y	
	Power Failure	Y	Y	Y	Y	
	High/Low Voltage	Y	Y	Y	Y	
	Sensor Error	Y	Y	Y	Y	
	Low Battery	Y	Y	Y	Y	
	High Ambient Temperature	Y	Y	Y	Y	
	Door Ajar	Y	Y	Y	Y	
	Start-up Delay	Y	Y	Y	Y	
	Accessories	Caster	Y	Y	Y	Y
Foot	Y	Y	Y	Y		
Porthole	Y/2	Y/2	Y/2	Y/2		
Shelves/Inner Doors	3/4	3/4	3/4	3/4		
USB Interface	Y	Y	Y	Y		
Remote Alarm	Y	Y	Y	Y		
5V Power Supply Port	Y	Y	Y	Y		
Temperature Chart Recorder	Optional	Optional	Optional	Optional		
RS485 Port	Y	Y	Y	Y		
CO <sub>2</sub> Backup System	Optional	Optional	Optional	Optional		
LN <sub>2</sub> Backup System	Optional	Optional	Optional	Optional		
Certifications	CE	Y	Y	/ Y	/ Y	
	UL	Y	Y	Y /	Y /	
	ENERGY STAR	Y	Y	Y Y	Y Y	

Haier Biomedical reserves the right to change products and specifications without prior notice.

Haier Biomedical's SmartFrequency technology manages the compressor speeds on demand. The ultra-low temperature freezers can thus achieve the world-leading energy efficiency in ultra-low temperature.



DW-86L829BP

## Smarter by Design

The SmartFrequency Conversion range of ultra low temperature freezers have been designed and developed at the Haier Biomedical R&D Institute. They are certified by one of the world's leading energy conservation and environmental protection organisations.



### Intelligent frequency Conversion technology

Two variable speed compressors are controlled for optimal freezer performance. Low energy consumption is achieved.



### Maximum energy efficiency

Our SmartFrequency Technology, coupled with our environmentally safe and friendly hydrocarbon refrigeration system, allows the Haier freezers to operate at a low level of energy of 8.2 kWh/day.



### Precise temperature control

The innovative control algorithm balances the effects of temperature loss with the unique frequency conversion refrigeration system, ensuring the cabinet temperature stability of  $\pm 3^{\circ}\text{C}$ .



### Frequency conversion Adaptive technology

Variable speed compressors in Haier Biomedical freezers are operated to produce the capacity that matches the demand of the load. The control system automatically tunes the speed of the compressors to optimize the operation.



### World-leading energy saving Refrigeration technology

Our hydrocarbon refrigeration technology can save energy by 50%, significantly reducing operator's cost. The refrigerants do not contain fluoride and chloride. The global warming potential is extremely low at 3. Thus they are very friendly to the environment.



### Minimal sound level output

Adaptive control technology controls the fans and compressors to reduce the sound level to 43.5 dB(A).



### Energy Saving Refrigeration

High efficiency cooling fans and compressors, combined with hydrocarbon refrigerants, ensure energy savings and long-term sample security



### Pressure Equalisation Port

Heated port with spring assisted mechanism to prevent icing on the vent allows users to reopen the main door soon after entering. Adopts chromium plating, rust-proof



### USB Interface

Enables users to download historical temperature data for compliance/auditing purposes



### Multilayered Sealing Structure

Triple layer of gaskets split between main and inner doors decreases heat loss and guarantee excellent warm up times in the event of a power failure



### Improved Handle Design

Lockable handle with unique key prevents other Haier freezer owner's access to your precious samples, also comes with space for a padlock for that extra security




### Multilevel Alarms

Alarming functions that include high and low temperature, sensor error, power failure, high ambient, clean filter and door ajar



DW-86L829BP

### DW-86L829BP | TYPICAL PERFORMANCE CHARACTERISTICS AT 25°C AMBIENT

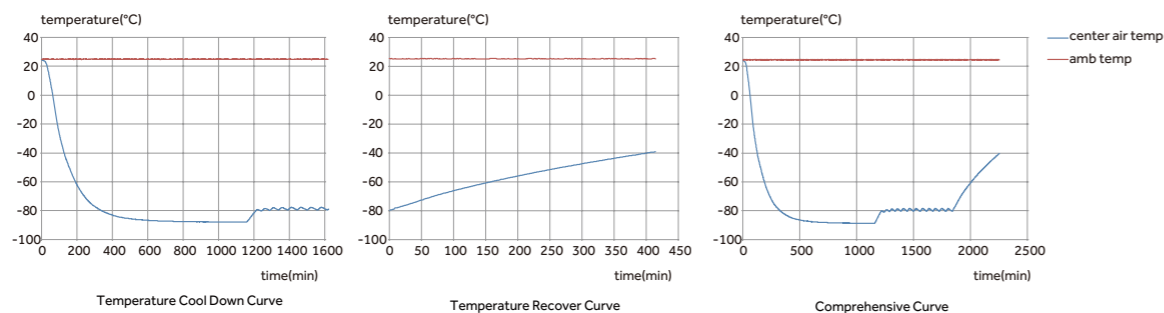

## Specifications

	Model	DW-86L579BP	DW-86L729BP	DW-86L829BP	DW-86L959BP	
Technical Data	Cabinet Type	Upright	Upright	Upright	Upright	
	Climate Class	N	N	N	N	
	Cooling Type	Direct Cooling	Direct Cooling	Direct Cooling	Direct Cooling	
	Defrost Mode	Manual	Manual	Manual	Manual	
	Refrigerant	HC	HC	HC	HC	
Performance	Sound Level (dB(A))	43.5	43.5	43.5	47	
	Cooling Performance (°C)	-86	-86	-86	-86	
Control	Temp Range (°C)	-40~-86	-40~-86	-40~-86	-40~-86	
	Controller	Microprocessor	Microprocessor	Microprocessor	Microprocessor	
Electrical Data	Display	LED	LED	LED	LED	
	Power Supply (V/Hz)	100~230/50/60	100~230/50/60	100-230/50/60 208-230/50/60	100-230/50/60 208-230/50/60	
Construction	Electrical Current (A)	14	14	15	6	
	Capacity (L/Cu.Ft)	579/20.4	729/25.7	829/29.2	959/33.9	
	Net/Gross Weight (approx)	kg	325/355	350/385	380/415	450/485
		lbs	716.5/782.6	771.6/848.8	837.7/914.9	992.1/1069.2
	Interior Dimensions (W*D*H)	mm	620*716*1310	766*716*1310	870*716*1310	1016*716*1310
		in	24.4*28.2*51.6	30.2*28.2*51.6	34.3*28.2*51.6	40.0*28.2*51.6
Exterior Dimensions (W*D*H)	mm	895*998*1980	1046*998*1980	1145*998*1980	1296*998*1980	
	in	35.2*39.3*78.0	41.2*39.3*78.0	45.1*39.3*78.0	51.0*39.3*78.0	
Packing Dimension (W*D*H)	mm	975*1095*2130	1120*1090*2130	1240*1080*2130	1410*1160*2185	
	in	38.4*43.1*83.9	44.1*42.9*83.9	48.8*42.5*83.9	55.5*45.7*86.0	
Loading Quantities	Container Load (20'/40'/40'H)	12/24/24	10/22/22	5/11/11	5/10/10	
	Remote Alarm (Dry contacts)	Y	Y	Y	Y	
Alarms	High/Low Temp	Y	Y	Y	Y	
	Hot Condenser	Y	Y	Y	Y	
	Power Failure	Y	Y	Y	Y	
	Sensor Error	Y	Y	Y	Y	
	Low Battery	Y	Y	Y	Y	
	High Ambient Temp	Y	Y	Y	Y	
	Door Ajar	Y	Y	Y	Y	
	Caster	Y	Y	Y	Y	
Accessories	Foot	Y	Y	Y	Y	
	Porthole	Y/2	Y/2	Y/2	Y/2	
	Shelves/ Inner Doors	3/4	3/4	3/4	3/4	
	USB Interface	Y	Y	Y	Y	
	5V Power Supply Port	Y	Y	Y	Y	
	Temp Recorder	Optional	Optional	Optional	Optional	
	RS485 Port	Y	Y	Y	Y	
	CO <sub>2</sub> Backup System	Optional	Optional	Optional	Optional	
	LN <sub>2</sub> Backup System	Optional	Optional	Optional	Optional	
	Certifications	CE	Y	Y	/	Y
UL		Y	Y	Y	/	
ENERGYSTAR		Y	Y	Y	Y	

BP suffix - Variable frequency inverter compressor  
Product appearance and specifications are subject to change without notice

Haier water-cooled ULT freezers are designed to meet strict requirements for storage of plasma, biological materials, vaccines, reagents, specimens, and other valuable samples. They are ideal for installations in hospitals, clinics, blood banks and medical research facilities where freezer heat must be removed by cooling water.



**Advantages**

- Intelligent frequency conversion technology
- Advanced control Touchscreen
- Higher efficiency to yield more energy savings of up to 20%
- Approximately 90% of the compressor heat generated during operation is removed by the cooling water, thus posing little impact to a laboratory's ambient
- Use less air-conditioning power for a comfortable laboratory condition
- Advanced control
- Low sound level
- Smart coolant control
- High quality industry grade hermetically sealed compressors
- Pressure protection system due to lack of water flow



**Specifications**

Model		DW-86L828W	DW-86L959W	
Technical Data	Cabinet Type	Upright	Upright	
	Climate Class	N	N	
	Cooling Type	Direct cooling	Direct cooling	
	Defrost Mode	Manual	Manual	
	Refrigerant	HC	HC	
	Sound Level (dB(A))	43.5	47	
Performance	Cooling Performance (°C)	-86	-86	
	Temperature Range (°C)	-40~-86	-40~-86	
Control	Controller	Microprocessor	Microprocessor	
	Display	LCD Touch Screen	LCD Touch Screen	
Electrical Data	Power Supply (V/Hz)	208-230/50/60	208-230/50/60	
	Electrical Current (A)	6	7	
Dimensions	Capacity (L/Cu.Ft)	828/29.2	959/33.9	
	Net/Gross Weight (approx)	kg	380/415	450/485
		lbs	837.7/914.9	992.1/1069.2
	Interior Dimension (W*D*H)	mm	870*716*1310	1016*716*1310
		in	34.3*28.2*51.6	40.0*28.2*51.6
	Exterior Dimension(W*D*H)	mm	1145*998*1980	1296*998*1980
		in	45.1*39.3*78.0	51.0*39.3*78.0
	Packing Dimension (W*D*H)	mm	1240*1080*2130	1410*1160*2185
		in	48.8*42.5*83.9	55.5*45.7*86.0
	Container Load (20'/40'/40'H)		5/11/11	5/10/10
Alarms	High/Low Temperature	Y	Y	
	Hot Condenser	Y	Y	
	Power Failure	Y	Y	
	High/Low Voltage	Y	Y	
	Sensor Error	Y	Y	
	Low Battery	Y	Y	
	High Ambient Temperature	Y	Y	
	Door Ajar	Y	Y	
	Caster	Y	Y	
Accessories	Foot	Y	Y	
	Porthole	Y/2	Y/2	
	Shelves/Inner doors	3/4	3/4	
	USB Interface	Y	Y	
	Remote Alarm (Dry contacts)	Y	Y	
	5V Power Supply Port	Y	Y	
	Temperature Recorder	Optional	Optional	
	Rs485 Port	Y	Y	
	CO <sub>2</sub> Backup System	Optional	Optional	
	LN <sub>2</sub> Backup System	Optional	Optional	
Others	Certification	CE	CE	

Product appearance and specifications are subject to change without notice


Haier Biomedical Salvum brand of ULT freezers has been designed to deliver energy savings and reduced carbon footprint. This range uses environmentally safe hydrocarbon refrigerants and high efficiency fan motors to maximize the cooling ability of the system and reduce energy consumption. While providing sample safety, the freezer design makes energy savings possible for laboratories.



DW-86L829



DW-86L579

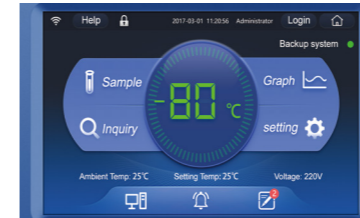


### Large Intelligent LCD Screen Visual Management


10-inch high-performance LCD capacitive screen, sensitive touch operation; intuitive display of inside temperature, ambient temperature, input voltage and other data and temperature curves.

### Inventory Management System

The single version of the inventory management allows users who do not have a racking system within the freezer to easily record item locations, entry and exit records.

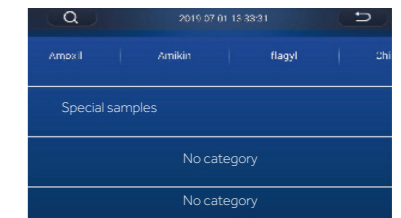

Single version (standard)

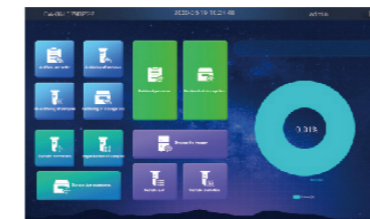



Professional version (optional)

In instances where customers use racking systems, the optional professional version of the inventory management is available. Freezer racks, boxes, and vials can easily record item locations and entry and exit records, facilitate item inventory and statistics, realize multi-screen interaction and reduce errors through secondary verification.

The freezer comes with the single version as standard.

### Equipped with Multiple Interfaces

- Standard USB interface, capable of storing data for 10 years
- Standard remote alarm terminal RS485 port

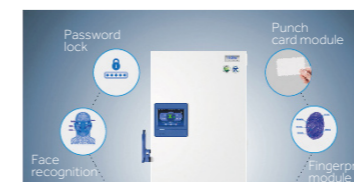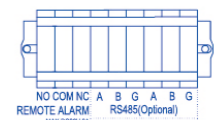

### Multiple Options for Safe Management and Control Systems

Standard password lock, optional fingerprint module, punch card module, face recognition, realizing safe and secure multi-user management.


### Optional IoT System Realizes Real-time Monitoring

The IoT APP can monitor the operating status of the equipment anytime and anywhere. Equipped with multiple alarm functions and able to self-diagnose system faults to ensure sample safety. (Contact your local area representative to check local availability.)



### Safe and Secure with Multiple Alarms

Multiple alarm functions include high temperature alarm, low temperature alarm, sensor failure alarm, power failure alarm, low battery alarm, door ajar alarm and high ambient temperature alarm.

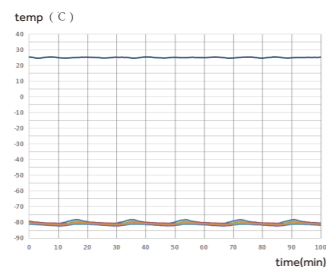

### Precise Temperature Control Ensures Safe Storage

Adjustable temperature set-point between -40°C ~ -86°C and the ability to monitor temperature curves in real time, equipped with sophisticated temperature safety alarm system to ensure product safety.

### Ergonomic Design



#### Narrow Door Design

The freezer can pass through an 750mm wide door when the door is open



#### New Handle Design

Adopting exclusive slide design, easy to open and close the door



#### Pressure Equalization Port

Equipped with pressure equalization port, realize easy opening  
Adopts chromium plating, rust-proof



#### Porthole

Two portholes as standard, allows for independent testing of cabinet temperature



#### Filter Design

Tool-free removable filter design, easy and quick to clean the filter



#### Detachable Inner Door Design

Tool-free detachable inner door, realizing easy and quick defrost

### Specifications

Model	DW-86L419	DW-86L579	DW-86L729	DW-86L829		
Technical Data	Cabinet Type	Upright	Upright	Upright	Upright	
	Climate Class	N	N	N	N	
	Cooling Type	Direct cooling	Direct cooling	Direct cooling	Direct cooling	
	Defrost Mode	Manual	Manual	Manual	Manual	
	Refrigerant	HC	HC	HC	HC	
Performance	Sound Level (dB(A))	47	47	52	52	
	Cooling Performance (°C)	-86	-86	-86	-86	
Control	Temp Range (°C)	-40~-86	-40~-86	-40~-86	-40~-86	
	Controller	Microprocessor	Microprocessor	Microprocessor	Microprocessor	
Electrical Data	Display	LCD Touchscreen	LCD Touchscreen	LCD Touchscreen	LCD Touchscreen	
	Power Supply (V/Hz)	220-240/50	220-240/50	220-240/50	220-240/50	
Construction	Electrical Current (A)	5	5.5	9	10	
	Capacity (L/Cu.Ft)	419/14.8	579/20.45	729/25.75	829/29.3	
Loading Quantities	Net/Gross Weight (approx)	kg	255/286	300/330	350/385	380/410
		lbs	562.2/630.5	661.4/727.5	771.6/848.8	837.7/903.9
	Interior Dimensions (W*D*H)	mm	465*716*1310	620*716*1310	766*716*1310	870*716*1310
		in	18.3*28.2*51.6	24.4*28.2*51.6	30.2*28.2*51.6	34.3*28.2*51.6
	Exterior Dimensions (W*D*H)	mm	830*980*1980	903*980*1960	1049*980*1980	1153*980*1980
		in	32.7*38.6*78.0	35.6*38.6*77.2	41.3*38.6*78.0	45.4*38.6*78.0
Packing Dimension (W*D*H)	mm	890*1090*2130	975*1095*2130	1120*1090*2130	1240*1080*2130	
	in	35.0*42.9*83.9	38.4*43.1*83.9	44.1*42.9*83.9	48.8*42.5*83.9	
Alarms	Container Load (20'/40'/40'H)	12/26/26	12/24/24	10/22/22	5/11/11	
	High/Low Temp	Y	Y	Y	Y	
	Hot Condenser	Y	Y	Y	Y	
	Power Failure	Y	Y	Y	Y	
	High/Low Voltage	Y	Y	Y	Y	
	Sensor Error	Y	Y	Y	Y	
	Low Battery	Y	Y	Y	Y	
	High Ambient Temp	Y	Y	Y	Y	
Door Ajar	Y	Y	Y	Y		
Accessories	Caster	Y	Y	Y	Y	
	Foot	Y	Y	Y	Y	
	Porthole	Y/2	Y/2	Y/2	Y/2	
	Shelves/ Inner Doors	3/2	3/4	3/4	3/4	
	USB Interface	Y	Y	Y	Y	
	Remote Alarm (Dry contacts)	Y	Y	Y	Y	
	5V Power Supply Port	Y	Y	Y	Y	
	Temp Recorder	Optional	Optional	Optional	Optional	
	RS485 Port	Y	Y	Y	Y	
	CO <sub>2</sub> Backup System	Optional	Optional	Optional	Optional	
LN <sub>2</sub> Backup System	Optional	Optional	Optional	Optional		
Certification	CE	Y	Y	Y	Y	
	Energy star	/	Y	Y	Y	

Product appearance and specifications are subject to change without notice

This product line is designed and manufactured for long term storage of various biological products, including viruses, germs, erythrocytes and leucocytes. Applications can be found in blood banks, hospitals, epidemic prevention services, research institutes, biological engineering institutes, laboratories in electronic and chemical plants.



**Advantages**

- World leading energy-efficient
- Hydrocarbon refrigeration system
- Slim cabinet design
- Reliable sample protection
- Malfunction alarms
- Excellent insulation performance



**Insulation and System Design**

- Special V-I-P (Vacuum Insulation Panel) insulation system reduces the heat gain by 25%
- High efficiency HC refrigeration system improves the overall efficiency by 45%
- Four individual insulated inner doors reduce the cold air loss to the minimum
- Heated Pressure Equalization Port makes re-accessing the unit fast
- About 50 dba sound level

**Safe and Reliable Storage**

- Superior temperature uniformity
- Dependable fans, compressors and other system related components




**Alarms (Visual and Audible)**

- Adjustable High/Low temperature alarm
- Sensor error
- Low battery
- Door ajar
- Power failure
- Hot condenser
- High ambient
- Remote alarm contact

### Extended Warm up Time During Power Failure

- Warm up time measures the time taken for temperature to rise up (from -80°C to -50°C) at 25°C ambient when the power is interrupted.
- Haier has the longest warm up time when compared with other major brands in the market.

TYPICAL PERFORMANCE CHARACTERISTICS AT 25 °C AMBIENT

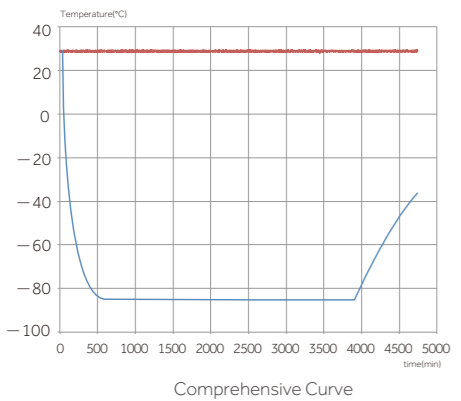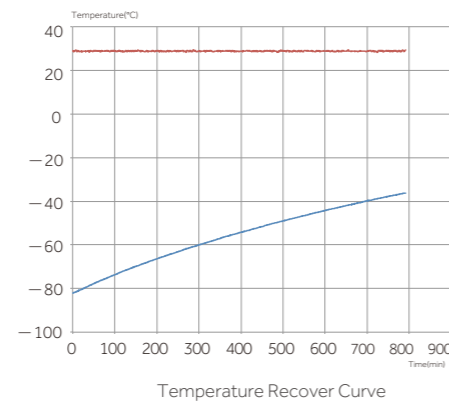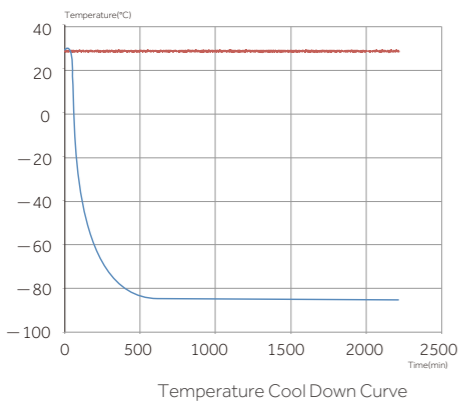

#### 4 Individual Removable Foam Inner Doors

- 4 individual inner doors can be opened independently to minimize frost buildup inside the chamber.
- Unique door seal design for the minimum loss of cold temperature during a door opening.
- Compatibility with existing racking system from competitors.
- Stainless steel handle to ensure proper strength for door latching.
- Some interior door handles have been upgraded to stainless steel for more comfortable door opening experience



#### Excellent Doors Seals

- Total of five gaskets to safeguard the freezer temperature, including four seals for the exterior door, one for each inner door



- Optimized slideway design, easier to open and close the door



#### Pressure Equalization Port

- Heated port with spring-assisted mechanism to prevent icing on the vent.
- Allows re-accessing the freezer after initial door opening.
- Adopts chromium plating, rust-proof



#### Circular- chart Recorder (Optional)


- Front-mounted.
- For independent temperature monitoring.



#### High Efficiency Refrigeration Components




Two port holes for ease of temperature monitoring

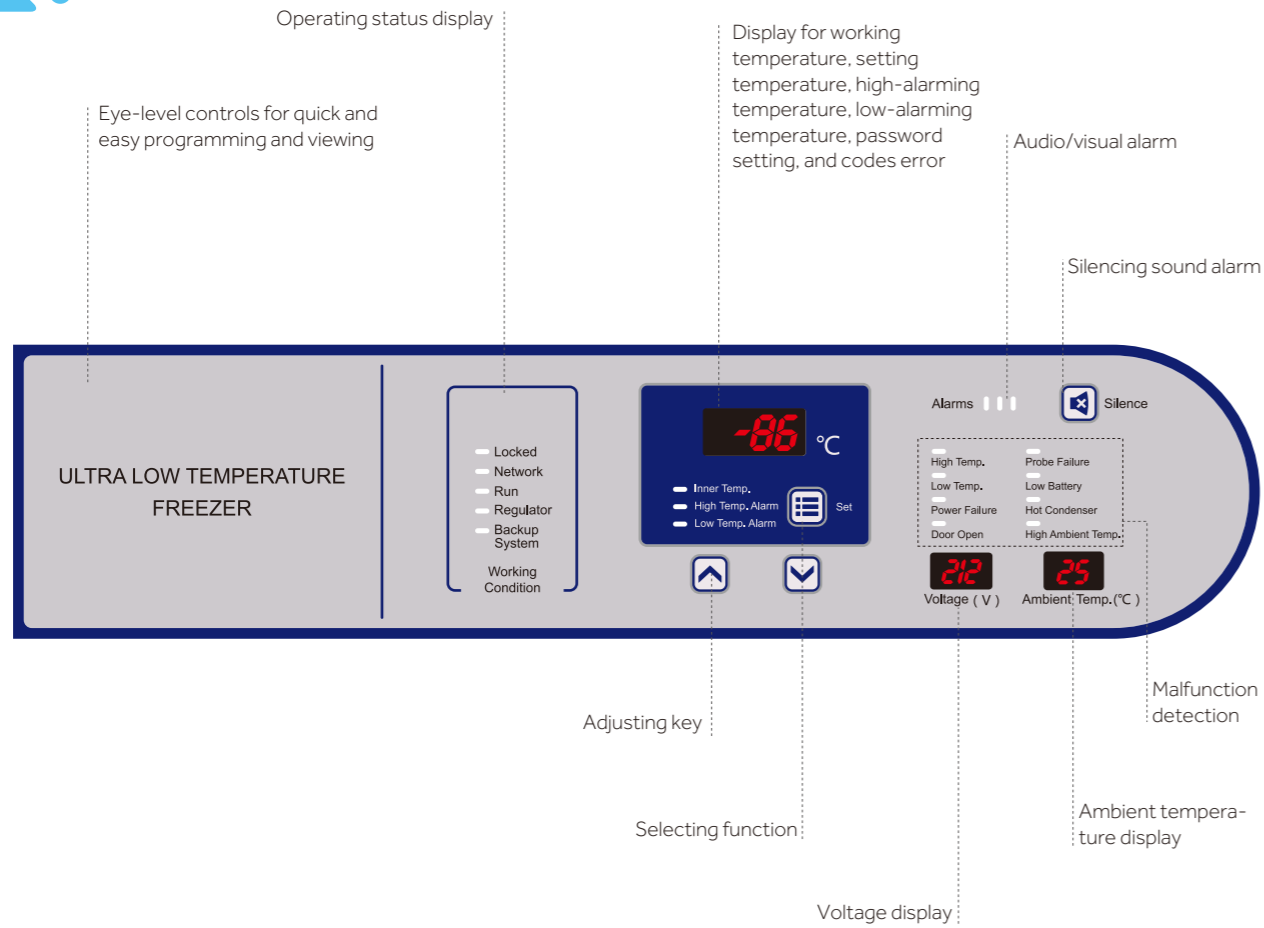


### Specifications

Alarm	Alarm Triggering Condition
High Temperature	Temperature reaches the warm alarm limit
Low Temperature	Temperature reaches the low alarm limit
Power Failure	Equipment loses power
Door Ajar	Door opening time exceeds set time, settable between 0 and 20 minutes
Sensor Error	E0.Ambient sensor fails E1.Condenser sensor fails E2.Main cabinet temperature control sensor fails E3.Heat exchanger sensor fails E4.Heat exchanger temperature fails
Low Battery	Battery capacity runs low or battery switch is not turned on
Hot Condenser	1. Condensers filter element is clogged 2. Ambient temperature is too high
High Ambient Temperature	Ambient temperature exceeds 32°C

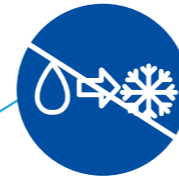

**Field proven reliability**

- Unique insulated inner door design for four separate storage compartments to minimize frost buildup inside the chamber
- Specialized control system design for a well-balanced operation of cascade refrigeration system
- Positive field proven reliability record

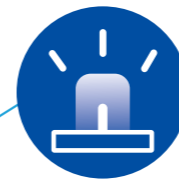

**Safety**

- Malfunction alarms including high and low temperature, power failure, sensor error, clean -filter, and extremely high ambient
- Capable of producing two types of alarm outputs: audible buzzer and visible flashing light
- Multiple built-in system protection features including user-settable protection code for controls, user settable delay to start, voltage compensation system, and protection against extreme high voltages
- Door open feature standard and USB port for temperature data downloading standard on upright models
- Remote alarm contacts



**Installation & Application**

- Wide range operating voltage system from 185V to 260 V designed to allow units installed in areas with poor voltage condition
- Suitable for 10°C to 32°C ambient temperature
- Input voltage and ambient temperature shown simultaneously for ease of monitoring environmental conditions
- Robust door latch designed for secure door closing
- Compact casters for ease of maneuvering



**Low sound level**

- Specialized refrigeration system design using whisper quiet fan and compressors
- Freezer chassis designed to absorb vibration and sound




**Energy saving**

- Unique door seal design for the minimum loss of cold temperature during a door opening
- High performance VIP insulation panels to minimize cabinet heat gain and to improve temperature stability
- Patented cabinet insulation system designed for optimal performance of cold storage temperature and minimal frost buildup
- Unique design of independent insulated inner door systems for independent access of storage space to provide the maximum protection of stored samples



**Key design features**

- Microprocessor-controlled system designed for controllable range of -40°C to -86°C for cabinet space with 1°C increment
- Large LED display for cabinet temperature, set temperature, ambient temperature, and input voltage
- Settable high temperature and low temperature alarms
- Automatic clean-filter alarm and sensor error alert
- Adjustable storage shelf height
- Optional temperature recorder, storage racks and storage boxes

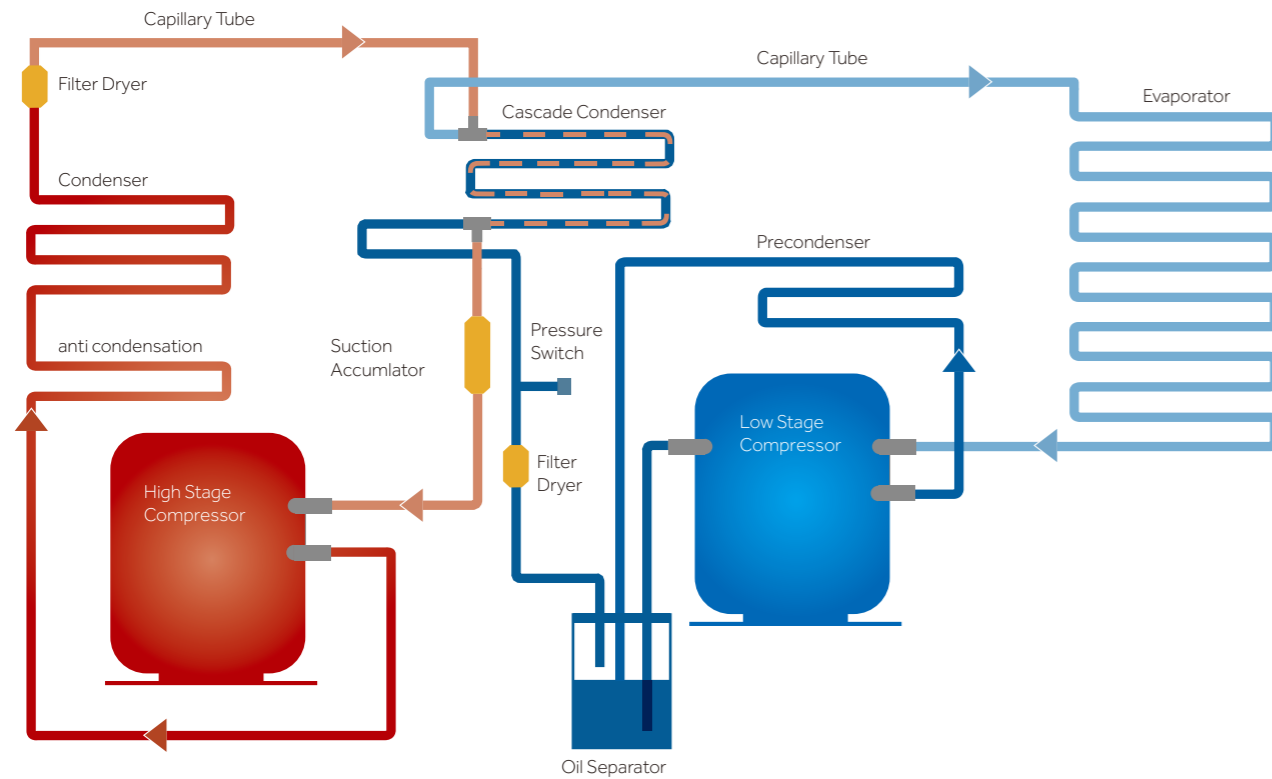

	<b>Filter dryer</b> • Filtration and drying device for water and debris
	<b>Capillary tube for high stage</b> • Expansion device for high stage liquid refrigerant
	<b>Capillary tube for low stage</b> • Expansion device for low stage refrigerant
	<b>High stage condenser</b> • Extra-large air cooled condenser dissipates product heat content efficiently.
	<b>Low stage evaporator</b> • Heat exchanger for expanded low stage refrigerant
	<b>Ultra low noise compressor</b> • Compressor provides reliable performance with ultra low noise • Specially designed low stage evaporator yields excellent temperature uniformity and recovery after loading and door opening. • Both high and low stage compressors are industrial grade hermetically sealed compressor designed for low temperature application. Sound level is extremely low.
	<b>Cascade condenser</b> • Heat exchanger between the high and low temperature stage in the system
	<b>Temprite oil separator</b> • Temprite oil separator can effectively separate oil and water to improve the refrigeration performance

Suitable for clinical, medical, scientific research, quarantine and other departments to store items under low temperature conditions. Applicable for universities, hospitals, disease prevention and control centers, blood stations, scientific research institutes, electronics and chemical enterprise laboratories and biomedical engineering research institutes. For storage of biological products and biological samples such as red and white blood cells, viruses, bone and bacteria. Also used for electronic devices and other materials used for cryogenic tests.



**DW-86L100J**



**Energy Efficient, Safe and Reliable**

High efficiency HC refrigeration system, optimised for energy efficiency delivering a power consumption figure of just 5.5kW/24hrs.



**Personal ULT Storage**

810mm cabinet height makes it easy to place on or under counter, saving storage space. Stackable design.



**Ergonomic design**

Ergonomic handle design ensures easy one-hand door opening.



**Low noise**

Optimized noise reduction cabinet and system design, emits sound level of only 46.8dB.

**VIP insulation and multilayered sealing design**

70mm insulating layer with 25mm VIP and 4 layers of gasket improves energy efficiency and reduces heat loss to deliver excellent warm up times in event of power failure.

**Optional IoT Module**

Real time monitoring of cabinet temperature, temperature setting, high and low temperature alarm value, temperature curve, alarm log and event log.

- User-settable parameters such as set point and alarms.
- Real-time cabinet temperature display, alarm information, power supply and compressor start/stop state.
- Standard USB port capable of storing >15 years of operating data for compliance.



Filter is easy to remove and clean without the need for tools.

4 casters + 2-foot locks, easy to move, lock and level.



Ergonomic design for easy door opening and closing. Lockable and equipped with 4 keys as standard with the ability to add a padlock for extra security when required.



Double stainless-steel inner doors to prevent cooling loss when opening the outer door, easy to clean.

**Microprocessor control system**

- Microcomputer electronic thermostat, LED temperature display, display precision 1°C, adjustable cabinet temperature set point -40°C~-86°C.
- Cabinet temperature/voltage/ambient temperature checking are available.
- Multiple alarm functions: high temperature alarm, low temperature alarm, sensor fault alarm, power failure alarm, low battery power alarm, open door alarm and high ambient temperature alert.
- Sound and light alarm mode, attachable to remote alarm interface.
- Battery backup alarm function operates continuously for >24hrs in the event of a power outage.
- Standard configuration: RS485 port and USB interface.
- Standard 5V power supply available for test equipment.
- Optional IoT module.

**Superior thermal insulation performance**

70mm super thick insulation layer design, aviation vacuum insulation material VIP, thickness of 25 mm or more, 4 layers of silicone seal, superior thermal insulation and energy saving effect.

**Porthole**

Portholes as standard, allows for independent testing of cabinet temperature.

**Security lock**

Standard door lock and padlock to ensure sample security and prevent unauthorised access.

**USB data storage**

Capable of storing more than 15 years of data.

**Specifications**

DW-86L486E



DW-86L490JA



DW-86L578JA



DW-86L628E

Model	DW-86L338J			DW-86L338JA			DW-86L388J			DW-86L486E			DW-86L490J		DW-86L490JA		DW-86L578J			DW-86L578JA			DW-86L628E																		
Technical Data	Cabinet Type	Upright									Upright						Upright																								
	Climate Class	N									N						N																								
	Cooling Type	Direct cooling									Direct cooling						Direct cooling																								
	Defrost Mode	Manual									Manual						Manual																								
	Refrigerant	HC									HC						HC																								
	Sound level (dB(A))	50									50						49																								
Performance	Cooling Performance (°C)	-86									-86						-86																								
	Temperature Range (°C)	-40~-86									-40~-86						-40~-86																								
Control	Controller	Microprocessor									Microprocessor						Microprocessor																								
	Display	LED									LED						LED																								
Electrical Data	Power Supply (V/Hz)	220~240/50			120/60			208~230/60			220~240/50			220~240/50			208~230/60		220~240/50			120/60			208~230/60			220~240/50													
	Electrical Current (A)	7.5			12			7.5			7.5			10			8		8		7.5			12			9			11											
Construction	Capacity (L/Cu.Ft)	338/11.9									388/13.7						486/17.1						490/17.3						578/20.4						626/22.1						
	Net/Gross Weight (approx)	kg	238/278									255/286						290/310						295/335						300/330						301/323					
		lbs	524.7/612.9									562.2/630.5						639.3/683.4						650.4/738.5						661.4/727.5						664.0/712.0					
	Interior Dimension (W*D*H)	mm	465*630*1165									465*716*1310						590*630*1310						590*630*1310						620*716*1310						760*630*1310					
		in	18.3*24.8*45.9									18.3*28.2*51.6						23.3*24.8*51.6						23.2*24.8*51.6						24.4*28.2*51.6						29.9*24.8*51.6					
	Exterior Dimension (W*D*H)	mm	823*900*1871									830*980*1980						953*900*1980						873*900*1980						903*980*1960						1035*900*1980					
in		32.4*35.4*73.6									32.7*38.6*78.0						37.5*35.4*78.0						34.4*35.4*78.0						35.6*38.6*77.2						40.7*35.4*78.0						
Packing Dimension (W*D*H)	mm	887*995*2040									890*1090*2130						1020*1020*2130						920*990*2130						970*1085*2130						1110*1005*2130						
	in	34.9*39.2*80.3									35.0*42.9*83.9						40.2*40.2*83.9						35.2*39.0*83.9						38.2*42.7*83.9						43.7*39.6*83.9						
Loading Quantities	Container load (20'/40'/40'H)	12/26/26									12/26/26						10/22/22						12/26/26						12/24/24						12/24/24						
Alarms	High/Low Temperature	Y									Y						Y						Y						Y												
	Hot Condenser	Y									Y						Y						Y						Y												
	Power Failure	Y									Y						Y						Y						Y												
	Sensor Error	Y									Y						Y						Y						Y												
	Low Battery	Y									Y						Y						Y						Y												
	High Ambient Temperature	Y									Y						Y						Y						Y												
	Door Ajar	Y									Y						Y						Y						Y												
Accessories	Caster	Y									Y						Y						Y						Y												
	Foot	Y									Y						Y						Y						Y												
	Porthole	Y/2									Y/2						Y/2						Y/2						Y/2												
	Shelves/Inner doors	3/2									3/2						3/4						3/4						3/4												
	USB Interface	Y									Y						Y						Y						Y												
	Remote Alarm (Dry contacts)	Y									Y						Y						Y						Y												
	5V Power Supply Port	Y									Y						Y						Y						Y												
	Temperature Recorder	Optional									Optional						Optional						Optional						Optional												
	RS232/485 Port	Optional									Optional						Optional						Optional						Optional												
CO <sub>2</sub> Backup System	Optional									Optional						Optional						Optional						Optional													
LN <sub>2</sub> Backup System	Optional									Optional						Optional						Optional						Optional													
Certifications	CE	Y			/			/			Y			/			/		Y			/			/			Y													
	UL	/			Y			Y			/			/			/		Y			/			Y			/													
	ENERGYSTAR	/			/			/			/			/			/		Y			Y			/			/													

**Specifications**

	Model	DW-86L728J			DW-86L728JA			DW-86L828J		DW-86L828JA		DW-86L100J		DW-86W100J		DW-86W420J		DW-86W420JA		
Technical Data	Cabinet Type	Upright						Upright						Upright		Chest		Chest		
	Climate Class	N						N						N		N		N		
	Cooling Type	Direct cooling						Direct cooling						Direct cooling		Direct cooling		Direct cooling		
	Defrost Mode	Manual						Manual						Manual		Manual		Manual		
	Refrigerant	HC						HC						HC		HC		HC		
	Sound level (dB(A))	50						50		51,5		46,8		49		50				
Performance	Cooling Performance (°C)	-86						-86						-86		-86		-86		
	Temperature Range (°C)	-40~-86						-40~-86						-40~-86		-40~-86		-40~-86		
Control	Controller	Microprocessor						Microprocessor						Microprocessor		Microprocessor		Microprocessor		
	Display	LED						LED						LED		LED		LED		
Electrical Data	Power Supply (V/Hz)	220~240/50		120/60		208~230/60		220~240/50		208~230/60		220~240/50		120/60		220~240/50		220~240/50   208~230/60		
	Electrical Current (A)	10		18		10		10		10		3		6,5		4		7,5		
Construction	Capacity (L/Cu.Ft)	728/25.7						828/29.2						100/3.5		100/3.5		420/14.8		
	Net/Gross Weight (approx)	kg	345/385		345/385		380/410		380/410		108/132		108/132		138/160		138/160		310/357	
		lbs	760.6/848.8		760.6/848.8		837.7/903.9		837.7/903.9		238/291		238/291		304.2/352.7		304.2/352.7		683.4/787.0	
	Interior Dimension (W*D*H)	mm	766*716*1310		766*716*1310		870*716*1310		870*716*1310		330*481*630		330*481*630		470*450*480		470*450*480		1367*462*652	
		in	30.2*28.2*51.6		30.2*28.2*51.6		34.3*28.2*51.6		34.3*28.2*51.6		13*19*25		13*19*25		18.5*17.7*18.9		18.5*17.7*18.9		53.8*18.2*25.7	
	Exterior Dimension (W*D*H)	mm	1049*980*1980		1049*980*1980		1145*980*1980		1145*980*1980		770*660*810		770*660*810		769*825*1120		769*825*1120		2130*870*1020	
in		41.3*38.6*78.0		41.3*38.6*78.0		45.1*38.6*78.0		45.1*38.6*78.0		30*26*32		30*26*32		30.3*32.5*44.1		30.3*32.5*44.1		83.9*34.3*40.2		
Packing Dimension (W*D*H)	mm	1120*1090*2130		1120*1090*2130		1240*1080*2130		1240*1080*2130		865*755*1020		865*755*1020		860*880*1310		860*880*1310		2230*930*1180		
	in	44.1*42.9*83.9		44.1*42.9*83.9		48.8*42.5*83.9		48.8*42.5*83.9		34.1*29.7*40.2		34.1*29.7*40.2		33.9*34.6*51.6		33.9*34.6*51.6		87.8*36.6*46.5		
Loading Quantities	Container load (20'/40'/40'H)	10/22/22		10/22/22		5/11/11		5/11/11		32/64/64		32/64/64		14/28/56		14/28/56		12/26/26		
Alarms	High/Low Temperature	Y						Y						Y		Y		Y		
	Hot Condenser	Y						Y						Y		Y		Y		
	Power Failure	Y						Y						Y		Y		Y		
	High/Low Voltage	/						/						/		/		/		
	Sensor Error	Y						Y						Y		Y		Y		
	Low Battery	Y						Y						Y		Y		Y		
	High Ambient Temperature	Y						Y						Y		Y		Y		
	Door Ajar	Y						Y						Y		Y		Y		
Accessories	Caster	Y						Y						Y		Y		Y		
	Foot	Y						Y						Y		Y		Y		
	Porthole	Y/2						Y/2						Y/1		Y/1		Y/1		
	Shelves/Inner doors	3/4						3/4						1/2		-/1		-/3		
	USB Interface	Y						Y						Y		Y		Y		
	Remote Alarm	Y						Y						Y		Y		Y		
	5V Power Supply Port	Y						Y						Y		N/A		N/A		
	Temperature Recorder	Optional						Optional						/		Optional		Optional		
	Rs232/485 Interface	Optional						Optional						NA/Y		Optional		Optional		
	CO <sub>2</sub> Backup System	Optional						Optional						Optional		Optional		Optional		
LN <sub>2</sub> Backup System	Optional						Optional						Optional		Optional		Optional			
Certifications	CE	Y		/		/		Y		/		Y		/		Y		/		
	UL	/		Y		Y		/		Y		/		Y		/		Y		
	ENERGY STAR	/		Y		/		/		Y		/		/		/		/		


Suitable for pharmaceutical, clinical and medical research use, for the long-term storage of samples and products under ultra-low temperature conditions.

**Advantage**



**Microprocessor Control System**

Microprocessor control, with temperature display precision at 0.1°C. Internal temperature is adjustable from of -40°C ~ -86°C

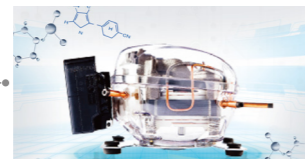

**Auto-cascade Hydrocarbon (HC) Refrigeration Technology**

Utilizing auto-cascade hydrocarbon compressor, with superior refrigeration and energy saving effect



**Excellent Doors Seals**

Total six gaskets to safeguard the freezer temperature, including a 70mm insulation layer, a 25mm vacuum VIP and 4 layers of silicone seal




**Stainless Steel Inner Door**

Prevents loss of internal air to keep samples safe



**Ergonomic Design**

Unique hook slide handle design, allows users to open and close the door with one hand

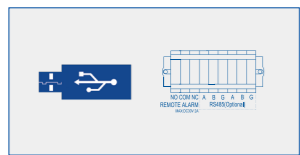

**Multiple Ports**

Enables users to download historical temperature data (up to 15 years) Equipped with remote alarm function and RS485 port



**Quick Access Filter Design**

Tool-free removable filter design, easy and quick to clean the filter



DW86L51J


**Internal Lock**  
Standard internal lock with 4 keys for multiple users



**Porthole**

One porthole as standard, allows for independent testing of cabinet temperature

**Advantage**

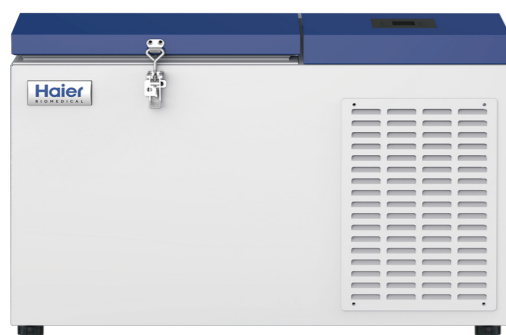

DW-80WZ15



**HC Energy Saving**

High quality HC high-efficiency compressor with optimized refrigeration system, which can save 50% on power consumption.



**One Unit for Multiple Applications**

\*Vehicle-mounted transportation with 24V(DC) power supply  
\*Long-term storage with 110V power supply.



**Superior Uniformity**

Built-in evaporator, fast cool down refrigeration, the temperature is safe, secure and reliable with the uniformity of ±3°C at key characteristic points.

**Flexible, Portable Mini-ULT Freezer**

**Ergonomic Design**



**Light Weight**

The whole unit weighs 27kg, easy to carry.



**Superior Insulation Performance**

25mm VIP + 70mm thick LBA foam insulation layer design, dual seal strips for better insulation.



**Safety Lock Catch**

Simple adjustable mechanical lock with lock catch, and it can be equipped with an external lock to guarantee the safety of sample transportation and storage.



**Metal Handle**

Hidden metal handle, stylish and comfortable, safe and durable, easy to carry.

**Specifications**

Model		DW-80WZ15		DW-86L51J	
Technical Data	Cabinet Type	Chest		Upright	
	Climate Class	N		N	
	Cooling Type	Direct Cooling		Direct Cooling	
	Defrost Mode	Manual		Manual	
	Refrigerant	HC		HC	
Performance	Sound Level (dB(A))	60		46.8	
	Cooling Performance (°C)	-80		-86	
	Temperature Range (°C)	-60~-80		-40~-86	
Control	Controller	Temperature Controller		Microprocessor	
	Display	LED		LED	
Electrical Data	Power Supply (V/Hz)	220/50, DC 24	110/60, DC 24	120/60	220~240/50
	Electrical Current (A)	2.5, 12	4.5, 12	7.5	3
Dimensions	Capacity (L/Cu.Ft)	15/0.53		51/1.8	
	Net/Gross Weight (approx)	(kg)	27/33	77/90	73/86
		(lbs)	59.52	169.8/198.4	160.9/189.6
	Exterior Dimension (W*D*H)	(mm)	690*320*500	532*640*806	
		(in)	27.2*12.6*19.7	20.9*25.2*31.7	
	Interior Dimension (W*D*H)	(mm)	260*180*330	330*481*316	
		(in)	10.2*7.1*13.0	13*16.5*12.4	
	Packing Dimension (W*D*H)	(mm)	805*430*725	570*735*1100	
		(in)	31.7*17.0*28.5	22.4*28.9*43.3	
	Cabinet Width (Without Handle and Hinge)	/		470	
Container Load (20'/40'/40'H)	105/210/210		60/126/126		
Alarms	High/Low Temperature	Y		Y	
	Hot Condenser	N		Y	
	Power Failure	Y		Y	
	High/Low Voltage	N		Y	
	Sensor Error	Y		Y	
	Low Battery	N		Y	
	High Ambient Temperature	N		Y	
	Door Ajar	N		Y	
	Start-up Delay	/		Y	
	Accessories	Caster	N		Y
Foot		Y		Y	
Porthole		N		Y/1	
Shelves/Basket		N/N		1/1	
USB Interface		N		Y	
Remote Alarm		N		Y	
Other	RS232/485 Port	N		NA/Standard	
	Certification	CE	UL	UL	CE

\*Haier Biomedical reserves the right to change products and specifications without prior notice.

Suitable for primary cells and cooling cells that are highly sensitive to temperature changes.

## Product Functions

**USB Interface**  
Easier to connect external modules such as printer and barcode reader

**Adjustable Feet**  
Convenient to adjust the height of machine

**Temperature Record**  
Information traceability: the frozen temperature curve and data can be transmitted to the computer for display, and PDF files can be generated synchronously for printing and archiving

**Handle**  
Easy to open and close the door, the user process is safer

**Lock**  
Ergonomic design, easy to operate

**Side-opening Door**  
Convenient access of frozen samples, suitable height

**Double Layer Sealing Structure**  
Good sealing performance

**Superior Insulation Performance**  
Environmental friendly polyurethane foam material has good thermal insulation performance

**Power Switch**  
Easy to operate, safe and reliable

**RS 232 Interface**

**Liquid nitrogen level monitoring port**

**Gas inlet**  
Pipe thread connection, safe and stable

**Pressure Relief Hole**  
When the inner pressure is too large, the pressure relief hole releases the pressure by itself to achieve the equilibrium state of internal and external pressure

**Power Connector**  
Plug and play

**Drain**


# Controlled Rate Freezer

## Product Advantages

- Precise Temperature Control**
  - Fine temperature control, the temperature rises and fall rate can be controlled in a large range, can achieve the cooling rate of 0.1~99°C/min
  - Stable temperature control, double solenoid valve design and special jet ring design, stable temperature control, effectively reduces cell damage
  - Precise temperature control, ADCR control algorithm, temperature fluctuation range is  $\pm 0.5^{\circ}\text{C}$
- Ergonomic Design**
  - 7-inch touch LCD screen, easy to operate
  - Liquid nitrogen allowance warning: Real-time display of liquid nitrogen allowance, timely warning when it is lower than 20%
  - Single handle press lock, simple and labor-saving, locking function, and the safety of the test process is guaranteed
  - Equipped with a real-time temperature curve printer and electronically archived temperature records

## Standard System Consists of

- Controlled Rate Freezer
- LN<sub>2</sub> Container
- LN<sub>2</sub> Transfer Hose
- Sample Probe
- Printer

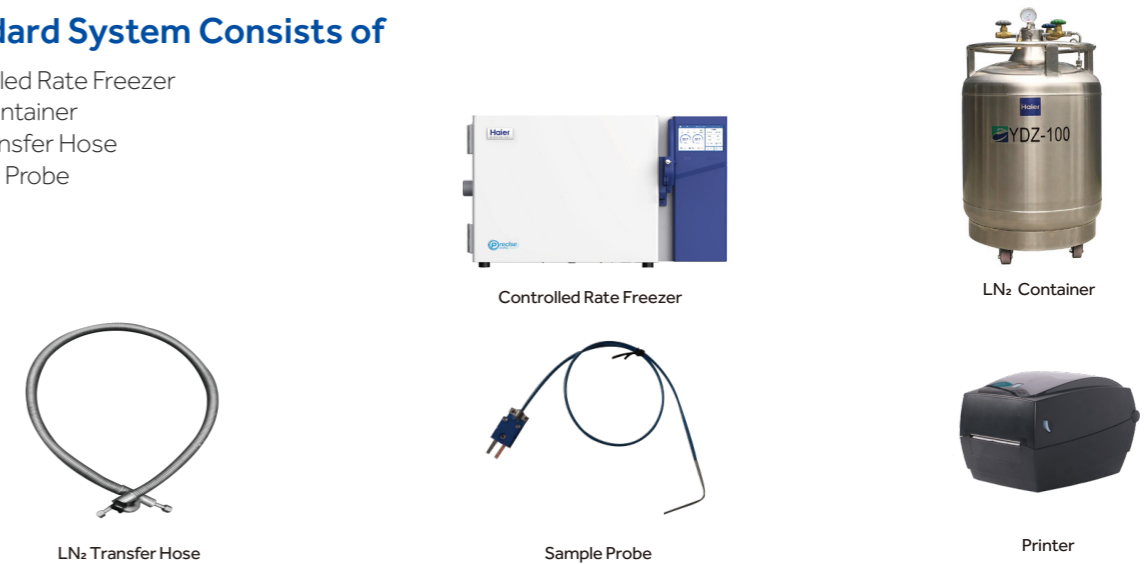

## Specifications

Model	Power Supply (V/Hz)	Rated Power (W)	Effective Volume (L)	Temperature Range (°C)
CJ-L37	230-50/60	1300	37	-180-50
Display Screen	Serial Port Form	Exterior Dimensions (W*D*H mm)	Interior Dimensions (W*D*H mm)	Net Weight (Kg)
7" colorful LCD touch screen	USB Interface/RS232 serial port	941*571*527	330/340/330	73

\*Haier Biomedical reserves the right to change products and specifications without prior notice.

### Advantage

- Small footprint, light and easy for installation
- User adjustable temperature setting
- Simple to program and operate
- Flexible – Can be installed on any ULT freezer, which has a port hole
- The CO<sub>2</sub> injection pipe is designed in combination of filter, to prevent the blockage of CO<sub>2</sub> back-up system

### Safety

- Liquid CO<sub>2</sub> test button to ensure the backup system is working
- Low CO<sub>2</sub> alarm system alerts the user when liquid CO<sub>2</sub> bottle is low in liquid level
- Ultra-low sound level compressors

### Reliability

- Stainless steel covering, more elegant
- Stainless steel input pipe design, allows for flexibility and ease of cylinder positioning
- Durable battery lasts up to 48 hours



### Specifications


	Model	HBX-IC
Technical Data	Climate Class	N
	Cooling Type	Direct cooling
	Refrigerant	CO <sub>2</sub>
Performance	Cooling Performance (°C)	-70
	Temperature Range (°C)	-40~-70
Control	Controller	Thermostat
	Display	LED
Electrical Data	Power Supply (V/Hz)	220-240/50
		208-230/60
	Power (W)	115/60
	Electrical Current (A)	20
Dimensions	Net/Gross Weight (approx)	11.2/14
		24.7/30.86
	Exterior Dimensions (W*D*H)	200*400*160
		7.8*15.7*6.3
	Packing Dimensions (W*D*H)	550*390*330
Alarms	CO <sub>2</sub> Margin Insufficient Alarm	21.7*15.4*13.0
	Low Battery Alarm	Y
	Sensor Error Alarm	Y
	Main Power Off Alarm	Y
	Charge Indicator	Y
	CO <sub>2</sub> Injection Test Button	Y
Accessories	Foot	4
Others	Certification	CE/UL

Product appearance and specifications are subject to change without notice



# LN<sub>2</sub> Backup System

LN<sub>2</sub> backup cooling system is an independent refrigeration system generally used for a -150°C mechanical cryo freezer. When there is a loss of power or the temperature of the freezer rises to the high alarm set point, the LN<sub>2</sub> backup system can be automatically activated to inject LN<sub>2</sub> into the chamber to maintain the freezer temperature. The backup system operates on a rechargeable battery when there is a loss of mains power.

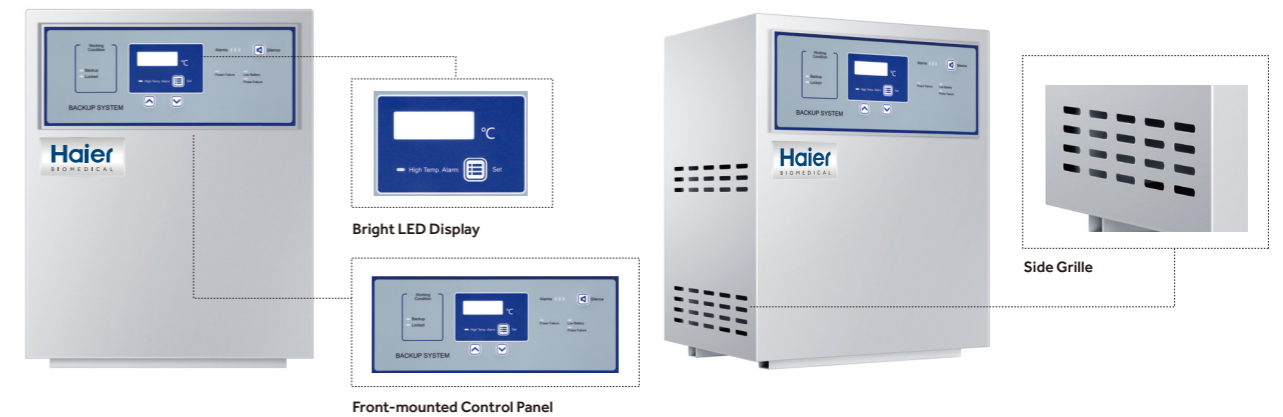

### Features

- Interlocked design to turn off the LN<sub>2</sub> injection when the door is opened
- Protection circuit to prevent over-charging battery
- Light weight and portable design, suitable for installation on top of a freezer
- Alarm functions to include low battery and sensor error

### Specifications

	Model	HBX-IIA
Technical Data	Cabinet Type	Horizontal
	Cooling Type	Direct cooling
	Refrigerant	LN <sub>2</sub>
Performance	Cooling performance (°C)	-135
	Temperature Range (°C)	-90~-135
Control	Controller	Thermostat
	Display	LED
Electrical Data	Power Supply (V/Hz)	220-240/50
		208-230/60
	Power (W)	115/60
Dimensions	Electrical Current (A)	20
		0.05
	Net/Gross Weight (approx)	15/17
		33.1/37.5
	Exterior Dimensions (W*D*H)	360*305*445
14.2*12.0*17.5		
Packing Dimensions (W*D*H)	450*375*720	
	17.7*14.8*28.3	
Alarms	Low Battery Alarm	Y
	Sensor Error Alarm	Y
	Main Power Off Alarm	Y
	Charge Indicator	Y
Others	Certification	CE/UL

Product appearance and specifications are subject to change without notice



**Consumables  
Solutions for  
ULT Freezer**

Haier Biomedical offers a full line of freezer storage racks to help making efficient use of space in chest and upright freezers, deep freezers, and ultra-low temperature freezers.

Normally, freezer racks may have horizontal access (for upright freezers) or vertical openings (for chest freezers). For optimal storage, client can choose racks designed to fit the particular freezer model or product.



## Product Features

- Manufactured using advanced capacitor energy storage resistance welding to product solid weld with smooth finish
- Various specifications for different storage requirements-frozen storage tube (box), SBS frozen storage tube (box), centrifugal tube, etc
- Handle is designed to reduce frost and ice build-up around handle
- Customised racking is available for different types of storage containers or models of freezer, contact us for more information

Rack Type	CAT.NO.	Material	Rack Layout (Rows*Columns)	Box Type	Maximum Box Dimensions(mm)	Rack Dimensions (D*W*H)(mm)	Net Weight (kg)
Sliding-drawer freezer rack	433001	SUS201	5*5	2 inch	133*133*53	686*141*291	5.1
	433002	SUS304	5*5	2 inch	133*133*53	686*141*291	5.1
	433003	SUS201	5*4	2 inch	133*133*53	552*141*291	4.3
	433004	SUS304	5*4	2 inch	133*133*53	552*141*291	4.3
	433005	SUS201	5*3	2 inch	133*133*53	417*141*291	3.8
	433006	SUS304	5*3	2 inch	133*133*53	417*141*291	3.8
	433007	SUS201	4*4	2 inch	133*133*53	552*141*240	3.1
	433008	SUS304	4*4	2 inch	133*133*53	552*141*240	3.1
	434040	SUS304	4*5	2 inch	133*133*53	686*141*233	缺参数
	434041	SUS304	3*3	2 inch	133*133*53	421*141*175	缺参数
	433009	SUS304	3*5	3 inch	133*133*94	686*141*291	3.1
433010	SUS304	3*4	3 inch	133*133*94	552*141*291	2.6	
433011	SUS304	7*5	For 0.5ml tube	133*133*37	686*141*291	6.8	
Side-access freezer rack	434036	SUS304	6*4	2 inch	133*133*53	558*140*320	6.3
	433012	SUS304	5*5	2 inch	133*133*53	686*141*291	5.1
	433013	SUS304	5*4	2 inch	133*133*53	552*141*291	4.3
	433037	SUS304	5*3	2 inch	133*133*53	417*141*291	3.92
	433014	SUS304	3*5	3 inch	133*133*94	686*141*291	3.1
	433015	SUS304	3*4	3 inch	133*133*94	552*141*291	2.7
433016	SUS304	4*4	2 inch	133*133*53	552*141*240	3.1	

# Freezer Storage Racks

Rack Type	CAT.NO.	Material	Rack Layout (Rows*Columns)	Box Type	Maximum Box Dimensions (mm)	Rack Dimensions (D*W*H)(mm)	Net Weight (kg)
SBS freezer rack	433017	SUS304	5*7	For 1.0ml/2ml sbs tube	86*128*50.8	686*141*291	6.8
	433018	SUS304	7*7	For 0.5ml sbs tube	86*128*38.5	686*141*291	6.2
	433019	SUS304	5*6	For 1.0ml/2ml sbs tube	86*128*50.8	552*141*291	5.2
	433020	SUS201	5*7	For 1.0ml/2ml sbs tube	86*128*50.8	686*141*291	5.3
	433021	SUS201	7*7	For 0.5ml sbs tube	86*128*38.5	686*141*291	6.2
	433022	SUS201	5*6	For 1.0ml/2ml sbs tube	86*128*50.8	552*141*291	5.2
Chest-type freezer rack	433023	SUS304	6 layer	2 inch	133*133*53	141*141*332	1.3
	433025	SUS304	8 layer	2 inch	133*133*53	141*141*442	1.4
	433026	SUS304	10 layer	2 inch	133*133*53	141*141*552	2
	433028	SUS304	12 layer	2 inch	133*133*53	141*141*662	2.3
Centrifugal tube freezer rack	433029	SUS304	5*22*2 layer	15ml centrifuge tube	/	686*141*291	3.66
	433030	SUS304	3*16*2 layer	50ml centrifuge tube	/	686*141*291	3.61

Haier Biomedical ULT Freezer Model	CAT.NO.	Rack and Box Discription	Freezer Rack Quantity	Freezer Box Quantity	Cryo Tube Quantity (2ml)
DW-86L338J/ DW-86L338JA	433013 (middle, 2 layers)	2 inch stainless steel side-access freezer rack 5*4 (2 inch freezer box*20 units)	6	120	12000
	433003/433004 (middle, 2 layers)	2 inch stainless steel sliding-drawer freezer rack 5*4 (2 inch freezer box*20 units)	6	120	12000
	433016 (top and bottom layers)	2 inch stainless steel side-access freezer rack 4*4 (2 inch freezer box*16 units)	6	96	9600
	433007/433008 (top and bottom layers)	2 inch stainless steel sliding-drawer freezer rack 4*4 (2 inch freezer box*16 units)	6	96	9600
DW-86L388J/ 418S/418ST/419	433012	2 inch stainless steel side-access freezer rack 5*5 (2 inch freezer box*25 units)	12	300	30000
	433001/433002	2 inch stainless steel sliding-drawer freezer rack 5*5 (2 inch freezer box*25 units)	12	300	30000
DW-86L486E	433013	2 inch stainless steel side-access freezer rack 5*4 (2 inch freezer box*20 units)	16	320	32000
	433003/433004	2 inch stainless steel sliding-drawer freezer rack 5*4 (2 inch freezer box*20 units)	16	320	32000
DW-86L490J/490JA	433013 (top 2 layers)	2 inch stainless steel side-access freezer rack 5*4 (2 inch freezer box*20 units)	8	160	16000
	433003/433004 (top 2 layers)	2 inch stainless steel sliding-drawer freezer rack 5*4 (2 inch freezer box*20 units)	8	160	16000
	433016 (bottom 2 layers)	2 inch stainless steel side-access freezer rack 4*4 (2 inch freezer box*16 units)	8	128	12800
	433007/433008 (bottom 2 layers)	2 inch stainless steel sliding-drawer freezer rack 4*4 (2 inch freezer box*16 units)	8	128	12800
DW-86L628E	433013	2 inch stainless steel side-access freezer rack 5*4 (2 inch freezer box*20 units)	20	400	40000
	433003/433004	2 inch stainless steel sliding-drawer freezer rack 5*4 (2 inch freezer box*20 units)	20	400	40000
DW-86L578J/578JA/ 578S/578ST/579/579BP/ 579BPT/578BPST	433012	2 inch stainless steel side-access freezer rack 5*5 (2 inch freezer box*25 units)	16	400	40000
	433001/433002	2 inch stainless steel sliding-drawer freezer rack 5*5 (2 inch freezer box*25 units)	16	400	40000
DW-86L728J/728JA/ 728S/728ST/728SAT/729/ 729BP/729BPT/728BPST	433012	2 inch stainless steel side-access freezer rack 5*5 (2 inch freezer box*25 units)	20	500	50000
	433001/433002	2 inch stainless steel sliding-drawer freezer rack 5*5 (2 inch freezer box*25 units)	20	500	50000
DW-86L828J/828JA/ 828S/828ST/829/829BP/ 829BPT/828W/828BPST	433012	2 inch stainless steel side-access freezer rack 5*5 (2 inch freezer box*25 units)	24	600	60000
	433001/433002	2 inch stainless steel sliding-drawer freezer rack 5*5 (2 inch freezer box*25 units)	24	600	60000
DW-86L959BP/ 959BPT/959W	433012	2 inch stainless steel side-access freezer rack 5*5 (2 inch freezer box*25 units)	28	700	70000
	433001/433002	2 inch stainless steel sliding-drawer freezer rack 5*5 (2 inch freezer box*25 units)	28	700	70000
DW-86L100J	433005/433006	2 inch stainless steel sliding-drawer freezer rack 5*3 (2 inch freezer box*15 units)	4	60	6000
DW-86W420J/420JA	433026	2 inch chest-type freezer rack (2 inch freezer box*10 units)	27	270	27000
DW-86W100J	433025	2 inch chest-type freezer rack (2 inch freezer box*8 units)	9	72	7200
DW-150W200	433026	2 inch chest-type freezer rack (2 inch freezer box*10 units)	12	120	12000
DW-150W209/ 150W209A	433026	2 inch chest-type freezer rack (2 inch freezer box*10 units)	12	120	12000
	433028	2 inch chest-type freezer rack (2 inch freezer box*12 units)	12	144	14400
DW-86L51J	433005/433006	2 inch stainless steel sliding-drawer freezer rack 5*3 (2 inch freezer box*15 units)	2	30	3000


Haier Biomedical can provide a range of high quality cryogenic storage tubes, which can withstand the extremes of cryogenic conditions while protecting the valuable samples. Meanwhile the Haier Biomedical cryogenic storage tubes have high compatibility with above freezer storage racks and ULT freezers, meeting the client's one-stop demands.



### Product Features

- Rigorous cytotoxicology-testing ensures pyrogen, RNase, DNase and endotoxins levels are less than 0.005EU/mL
- Complying with ANSI/SLAS international standards, suitable for a variety of automated storage and processing equipment
- Laser etching process and three-code information integration technology are adopted to ensure the security and redundancy to provide an extra layer of protection
- The snap design of the base is easy to use and prevents scratches after freezing
- Designed and manufactured to prevent leakage at -98kpa to meet the IATA dangerous goods PI650 standard

CAT.NO.	Description	Unit Quantity	Packing Specification
413001	0.5ml (ultra-low temperature volume 0.5ml), slender, external spin tube, bottom prefabricated QR code, side wall prefabricated digital code and 1D barcode; box side prefabricated digital code and 1D barcode, box bottom prefabricated QR code, 8*12 standard SBS layout	7680	96pcs/box, 80boxes/case
413002	0.75ml (ultra-low temperature volume 0.75ml), slender, external spin tube, bottom prefabricated QR code, side wall prefabricated digital code and 1D barcode; box side prefabricated digital code and 1D barcode, box bottom prefabricated QR code, 8*12 standard SBS layout	6720	96pcs/box, 70boxes/case
413003	1.0ml (ultra-low temperature volume 1.0ml), slender, external spin tube, bottom prefabricated QR code, side wall prefabricated digital code and 1D barcode; box side prefabricated digital code and 1D barcode, box bottom prefabricated QR code, 8*12 standard SBS layout	4800	96pcs/box, 50boxes/case
413004	2.0ml (ultra-low temperature volume 1.9ml), slender, external spin tube, bottom prefabricated QR code, side wall prefabricated digital code and 1D barcode; box side prefabricated digital code and 1D barcode, box bottom prefabricated QR code, 6*8 standard SBS layout	2400	48pcs/box, 50boxes/case
413005	2.0ml (ultra-low temperature volume 1.9ml), slender, external spin tube, bottom prefabricated QR code, side wall prefabricated digital code and 1D barcode; box side prefabricated digital code and 1D barcode, box bottom prefabricated QR code, 10*10 square layout	3000	100pcs/box, 30boxes/case
413006	0.5ml (ultra-low temperature volume 0.5ml), slender, external spin tube, bottom prefabricated QR code, side wall prefabricated digital code and 1D barcode	13200	220pcs/bag, 60bags/box
413007	1.0ml (ultra-low temperature volume 1.0ml), slender, external spin tube, bottom prefabricated QR code, side wall prefabricated digital code and 1D barcode	6000	150pcs/bag, 40bags/case
413008	2.0ml (ultra-low temperature volume 1.9ml), slender, external spin tube, bottom prefabricated QR code, side wall prefabricated digital code and 1D barcode	3500	70ticks/bag, 50bags/box
413009	Cryopreservation box, matching 0.5ml slender external rotation cryopreservation tube, prefabricated digital code and 1D barcode on the side of the box, prefabricated QR code at the bottom of the box, 8*12 standard SBS layout	80	80boxes/case
413010	Cryopreservation box, matching 0.75ml slender external rotation cryopreservation tube, prefabricated digital code and 1D barcode on the side of the box, prefabricated QR code at the bottom of the box, 8*12 standard SBS layout	70	70boxes/case
413011	Cryopreservation box, matching 1.0ml slender external rotation cryopreservation tube, prefabricated digital code and 1D barcode on the side of the box, prefabricated QR code at the bottom of the box, 8*12 standard SBS layout	50	50boxes/case
413012	Cryopreservation box, matching 2.0ml slender external rotation cryopreservation tube, prefabricated digital code and 1D barcode on the side of the box, prefabricated QR code at the bottom of the box, 6*8 standard SBS layout	50	50boxes/case
413013	Cryopreservation box, matching 2.0ml slender external rotation cryopreservation tube, prefabricated digital code and 1D barcode on the side of the box, prefabricated QR code at the bottom of the box, 10*10 standard SBS layout	30	30boxes/case

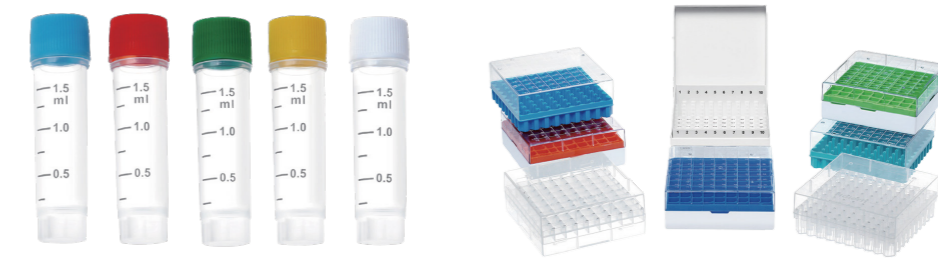

### Product Features


- Detachable color label design, color classification of samples according to requirements
- Clear scale and large writing area
- Internal rotation and external rotation are designed to meet different sample storage requirements
- Multi-color tube cap design, according to the needs of the sample color classification
- 2D tube bottom prefabricated two-dimensional code, laser etching technology, wear and corrosion resistant
- Pipe wall prefabricated bar code and digital code, three codes in one, to meet different scene identification

CAT.NO.	Description	Unit Quantity	Packing Specification
423001	Paper freezer box, 2 inches 81 compartments, white, flip lid	50	50boxes/case
423002	Paper freezer box, 2 inches 100 compartments, white, flip lid	50	50boxes/case
421003	81-well PP freezer box, -121-90°C, monochrome box, white, yellow, red, blue, green optional, flat lid	50	50boxes/case
421004	100-well PP freezer box, -121-90°C, monochrome box, white, yellow, red, blue, green optional, flat lid	50	50boxes/case
422005	25-well PC freezer box, -196-121°C, monochrome box, purple, red, blue, green optional, flat lid	50	50boxes/case
422006	81-well PC freezer box, -196-121°C, monochrome box, purple, red, blue, green optional, flat lid	50	50boxes/case
422007	100-well PC freezer box, -196-121°C, monochrome box, purple, red, blue, green optional, flat lid	50	50boxes/case
422010	3-inch 81-well PC freeze storage box, -196-121°C, single color one box, white, red, blue, green optional, flat lid	20	20boxes/case
422008	Supporting 0.5ml external rotation cryopreservation tube, 100-well PC cryopreservation box, -196-121°C, blue (supporting H-ES05B; H-ES05D-TS)	50	50boxes/case
422009	Supporting 1.5ml external rotation cryopreservation tube, 100-well PC cryopreservation box, -196-121°C, blue (supporting H-ES15B; H-ES15D-TS)	50	50boxes/case
423010	0.5ml cryopreservation tube, external rotation, self-standing bottom, with scale, with white writing area, -196°C - 121°C	1000	100pcs/bag, 10bags/case
423011	1.5ml cryopreservation tube, external rotation, self-supporting bottom, with scale, with white writing area, -196°C - 121°C	1000	100pcs/bag, 10bags/case
423012	2ml cryopreservation tube, external rotation, self-supporting bottom, with O-ring, with scale, with white writing area, -196°C - 121°C	1000	100pcs/bag, 10bags/case
423013	5ml cryopreservation tube, external rotation, self-supporting bottom, with O-ring, with scale, with white writing area, -196°C - 121°C	1000	100pcs/bag, 10bags/case
423014	1ml cryopreservation tube, internal rotation, self-supporting bottom, with O-ring, with scale, with white writing area, -196°C - 121°C	1000	100pcs/bag, 10bags/case
423015	2ml cryopreservation tube, internal rotation, self-supporting bottom, with O-ring, with scale, with white writing area, -196°C - 121°C	1000	100pcs/bag, 10bags/case
423016	0.5ml cryopreservation tube (boxed), external rotation, self-supporting bottom, with scale, with white writing area, -196°C - 121°C	1200	100pcs/box, 12boxes/case
423017	1.5ml cryopreservation tube (boxed), external rotation, self-supporting bottom, with scale, with white writing area, -196°C - 121°C	1200	100pcs/box, 12boxes/case
423018	0.5ml 2D cryopreservation tube (bagged), external rotation, self-supporting bottom, with scale, bottom QR code, side barcode + digital code, -196°C-121°C	1000	100pcs/bag, 10bags/case
423019	1.5ml 2D cryopreservation tube (bagged), external rotation, self-supporting bottom, with scale, bottom QR code, side barcode + digital code, -196°C-121°C	1000	100pcs/bag, 10bags/case
423020	0.5ml 2D cryopreservation tube (boxed), external rotation, self-supporting bottom, with scale, bottom QR code, side barcode + digital code, -196°C-121°C	1200	100pcs/box, 12boxes/case
423021	1.5ml 2D cryopreservation tube (boxed), external rotation, self-supporting bottom, with scale, bottom QR code, side barcode + digital code, -196°C-121°C	1200	100pcs/box, 12boxes/case


The Hoierh high speed batch scanner can be used to scan various coded, boxed prefabricated code cryogenic tubes to achieve high-speed identification and efficient management of samples.




Intelligent induction of tube boxes and automatic light filling function, adaptable to multiple applications

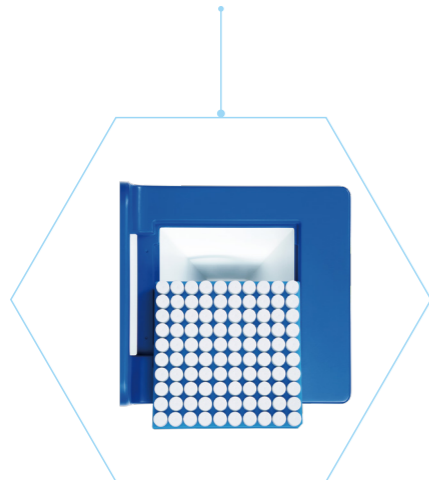

Reading synchronously the bottom code of freezing tubes and of boxes

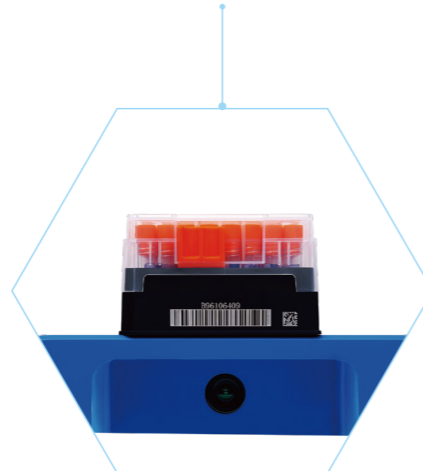

Detachable module design, compatible with SBS, square and other cryogenic boxes

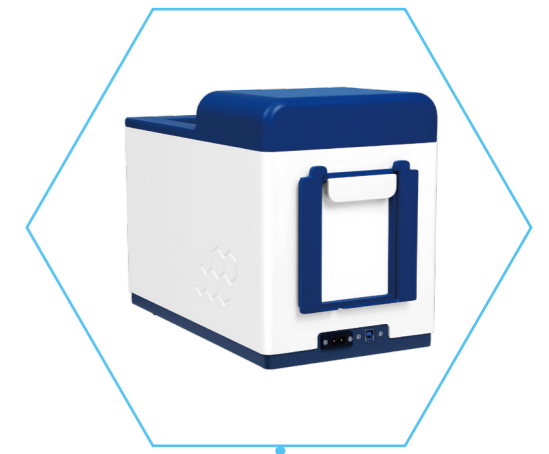

Ergonomic design: The scanner can be wall-mounted and stored safely using the attachment included

Classification	Model	Parameter	Sets/Case	Package
High speed batch scanner	HS-SBS-001	Batch scanning, high-speed scanning, side code scanning, suitable for 10*10/SBS and other specifications of prefabricated coded standard cryogenic tube/box of various brands	1	1 set/case
	HS-SBS-002	Batch scanning, high-speed scanning, side code scanning, suitable for 10*10/SBS and other specifications of prefabricated coded standard cryogenic tube/box of various brands	1	1 set/case
Top batch scanner	HS-100-002	Batch scanning, high-speed scanning, top scanning, suitable for scanning of 10*10/9*9 cryogenic tube, centrifuge tube, etc. (two-dimensional code label is pasted on the top)	1	1 set/case

# Single Tube Scanner

The Hoierh single tube scanner is flexible and quick to scan information



Wireless communication distance up to 30 meters



Support battery charging and replacement, can work for 20 hours or standby for 50 hours when fully charged



Excellent DPM code reading ability, suitable for all kinds of precision code recognition, fast and accurate

Classification	Model	Parameter	Sets/ case	package
Single tube scanner	HS-H-001	Handheld scanner, provide single cryogenic tube or cryogenic box, scanner recognition one-dimensional code/two-dimensional code function, plug and play	1	1 set/case
	HS-W-001	Wireless handheld scanner, provide single cryogenic tube or cryogenic box, scanner recognition one-dimensional code/ two-dimensional code function, plug and play	1	1 set/case
	HS-T-001	Desktop scanner, provide single cryogenic tube or cryogenic box, scanner recognition one-dimensional code/two-dimensional code function, plug and play	1	1 set/case

# Batch Cover Opener

The Hoierh batch cover opener opens or closes the multi-channel cap which provides increased efficiency compared to manual operation and reduces the risk of contact with biochemical reagents and other safety risks



Simultaneously open up to 8 channels at one time



Convenient and fast removal of cover



Support design, convenient storage opened batch tube cover

Classification	Model	Parameter	Sets/ case	package
Batch cover opener	BQ-OP8	The 8-channel hand-held electric cover opener can open the prefabricated QR code cryogenic tube (8x12) with 8 cryogenic tube covers in a row, or (diamond layout 8x6) with 4 cryogenic tube covers in a row	1	1 set/case
	BQ-OP4	The 4-channel hand-held electric cover opener can open the prefabricated QR code cryogenic tube (4x12) with 4 cryogenic tube covers in a row, or (4x6) with 4 cryogenic tube covers in a row	1	1 set/case

The low-temperature printer allows users to print different labels that are durable, stable, and resistant to wear, ensuring the long-term safety of sample information.



Suitable for all kinds of experimental instruments, resistant to solvent and chemical corrosion, repeated freezing and thawing without falling off

Touch operation for quick adjustment of printer Settings



High performance printing, meet the high requirements of biobank processing

Classification	Product Number	Model	Supporting Printer	Parameter	Loading Capacity
Printer	i5100 300dpi	Printer i5100 300dpi	\	Thermal transfer printer, printing accuracy 300dpi, maximum printing speed 300mm/s	1 set/case
	BBP12	Brady BBP12 printer	\	Printing accuracy 300dpi, low-end printing	1 set/case
	i5100 600dpi	Printer i5100 600dpi	\	Thermal transfer printer, printing accuracy 600dpi, maximum printing speed 150mm/s	1 set/case
Ribbon	IP-R4302	Ribbon IP R4302	i5100 300dpi	60mm*300m/roll, used together with IP 300; Print about 4-5 rolls of labels (B490A or B490B)	1 roll/bag
	IP-R4902	Ribbon IP-R4902	i5100 300dpi	60mm*300m/roll, used together with IP 300; Print about 4-5 rolls of labels (B6422)	1 roll/bag
	R4313	Ribbon R4313	BBP11/BBP12	110mm*74m/roll, used with BBP11/BBP12 printer	1 roll/bag
Label	THT-152-053-3-2A	B053 material (25.4*9.53mm)*2 two rows of labels	i5100 300dpi	Label for room temperature use. The label material is synthetic paper, and the use temperature range is -20°C~+80°C; Two labels in a row, 3000 pcs/roll, and the matching ribbon is IP-R4302	3000 pcs/roll
	THT-152-053-3-RP	B053 material 25.4*9.53mm one in a row	i5100 300dpi	Label for normal temperature use. The label material is synthetic paper, and the use temperature range is -20°C~+80°C; Two labels in a row, 4000pcs/roll, and the matching ribbon is IP-R4302; Customized, pre printed coding	4000 pcs/roll
	THT-163-6422-3	B6422 lable 25.4mm*9.53mm + D=9.53mm	i5100 300dpi	Apply to cryogenic tube below 2ml; Environment temperature above 40 °C	3000 pcs/roll
	THT-163-490-3	B490 material 25.4mm*9.53mm + D=9.53mm	i5100 300dpi	Apply to cryogenic tube below 2ml and pasted on the wall and round cover of the cryogenic tube; Can withstand can be directly pasted after taking out the cryogenic tube in liquid nitrogen at room temperature	3000 pcs/roll
	THT-152-490-3	B490 material 25.4mm*9.53mm	i5100 300dpi	Apply to cryogenic tube below 2ml, only pasted on the wall of cryogenic tube; Can withstand can be directly pasted after taking out the cryogenic tube in liquid nitrogen at room temperature	3000 pcs/roll
	THT-151-490-3	B490 material 15.24mm*10.16mm+D=8.26mm	i5100 300dpi	Apply to cryogenic tube below 2ml and pasted on the wall and round cover of the cryogenic tube; Can withstand can be directly pasted after taking out the cryogenic tube in liquid nitrogen at room temperature	3000 pcs/roll

Cryo Freezer



DW-150W209

-86°C ULT Freezer



DW-86L578S DW-86L578S(A)T DW-86L728S DW-86L728S(A)T  
 DW-86L579BP/729BP DW-86L829BP/959BP  
 DW-86L578BPST DW-86L728BPST DW-86L828BPST

Water Cooling Freezer



DW-86L579BPT/729BPT DW-86L829BPT/959BPT  
 DW-86L828W DW-86L959W  
 DW-86L419/579 DW-86L729/829 DW-86L338J  
 DW-86L338J DW-86L388J  
 DW-86L490J DW-86L490JA  
 DW-86L578J DW-86L578JA



DW-86L728J DW-86L728JA  
 DW-86L828J DW-86L828JA  
 DW-86L486E DW-86L628E  
 DW-86L51J DW-86L100J

Mini ULT Freezer

Chest




DW-86W100J DW-86W420J DW-86W420JA

Controlled Rate Freezer



DW-80WZ15

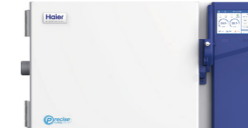

CJ-L37

CO<sub>2</sub> Backup System



HBX-IC HBX-IIA

Freezer Storage Racks



Barcode Cryogenic Tubes and Storage Boxes



Universal Cryogenic Tubes and Storage Boxes



High Speed Batch Scanner



Single Tube Scanner



Batch Cover Opener



Low Temperature Label Printer

